


PressureGuard® CFT®

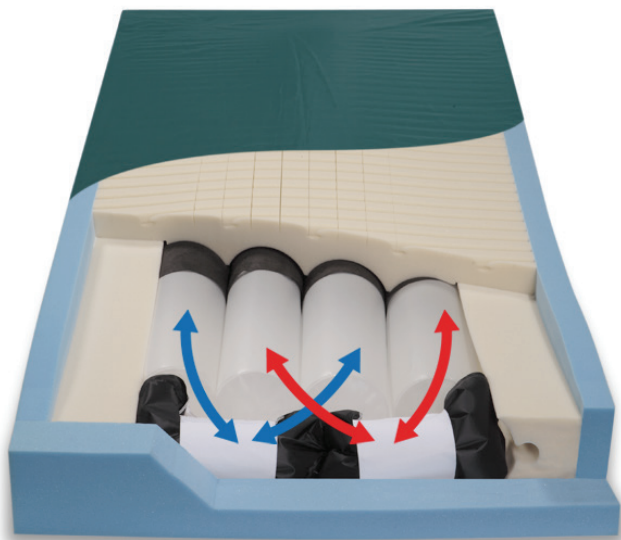


The revolutionary PressureGuard CFT's Constant Force Technology® automatically adjusts a network of interconnected air cylinders and elasticized reservoirs to the appropriate, therapeutic level, regardless of the user's weight or position.

The CFT non-powered system is approved for reimbursement as a Group II support surface, based on head-to-head studies in which it outperformed leading powered systems in ulcer healing.

The CFT is extremely cost effective. Requires virtually no set-up, programming, or maintenance other than periodic air boost (Air System Renewal Kit is required). Wipes clean. Generates no noise or heat and uses no electricity, making it immune to power interruptions.

Use it in-facility to slash rental costs, increase user satisfaction, and improve regulatory compliance on issues from skin care to bed safety. Or, select it as a trouble-free but supremely effective option for at-home users.

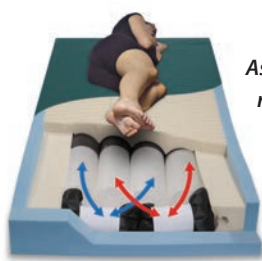

As user shifts position, reservoirs absorb displaced air to accommodate load and re-equalize pressures.

Therapy Provided

- ✓ Non-powered (Air) Reactive Pressure Redistribution

QUICK FACTS

Indications:

- Treatment of stage 1-4 pressure ulcers
- Prevention of skin breakdown in the high risk or severely compromised patient
- Patient can be positioned off the ulcer in at least 2 positions¹
- History of previously healed pressure ulcers
- Can be repositioned frequently
- Poorly padded bony prominences

Key Advantages:

- Proven for both ulcer healing and prevention
- Outperforms powered systems for wound healing*
- Medicare covered (non-powered Group II)
- No controls to adjust
- No reliance on power source
- PressureGuard® design: stability, safety, anti-shearing

Item #	Description
8075-29	PressureGuard CFT®, 75"L x 35"W x 7"H
* 8280-29	PressureGuard CFT®, 80"L x 32"W x 7"H
8080-29	PressureGuard CFT®, 80"L x 35"W x 7"H
8084-29	PressureGuard CFT®, 84"L x 35"W x 7"H
* CF8042-29	PressureGuard CFT®, 80"L x 42"W x 7"H Bariatric
* CF8048-29	PressureGuard CFT®, 80"L x 48"W x 7"H Bariatric
* CF8054-29	PressureGuard CFT®, 80"L x 54"W x 7"H Bariatric
* CF8442-29	PressureGuard CFT®, 84"L x 42"W x 7"H Bariatric
* CF8448-29	PressureGuard CFT®, 84"L x 48"W x 7"H Bariatric
* CF8454-29	PressureGuard CFT®, 84"L x 54"W x 7"H Bariatric
SERVICE-00	Air System Renewal Kit (includes hand pump and auto-adjust valve)

Twin/Full/Queen/King sizes available: See In-Home Styles – page 19.
Additional bariatric options available – see page 18.

* Allow up to 10 working days lead time for these sizes.

Medicare Approved: E0373 (Group II)
Weight limit: 500 lbs.
750 lbs. Bariatric models and Full/Queen/King sizes
Warranty: 5 years, non-prorated
Mattress Weight: 21 lbs.

All models are fire barrier version (FBV) as indicated by "-29" suffix, denoting certified conformance to NFPA 101 (Life Safety Code); Cal Tech Bulletin #117; 16 CFR Parts 1632 and 1633; and BFD IX-11.

¹In situations where positioning on the pressure ulcer cannot be avoided (e.g. the individual has multiple ulcers on multiple surfaces), limit the amount of time the individual is positioned on the ulcer. Source: NPAUP/EPJAP Pressure Ulcer Prevention and Treatment, Clinical Practice Guideline, p 66-67

*Studies available upon request

