

RAPID DIAGNOSTICS

osom[®]

Trichomonas Rapid Test

FOR THE DETECTION OF *TRICHOMONAS VAGINALIS* ANTIGENS

(May be used with patient collected vaginal swabs)



www.osomtests.com



FEATURES:

- ✓ **ACCURATE**; 83% sensitivity vs. culture; 95% agreement with culture and wet mount combined
- ✓ **FAST**; Results in 10 minutes or less
- ✓ **FLEXIBLE**; Contains patient instructions to allow self-collection of sample
- ✓ **OBJECTIVE**; Detects the antigen; removes uncertainty associated with live organism instability
- ✓ **TRAINING**; Online training and certification; www.osomtraining.com
- ✓ **CUSTOM CARE**; Technical support staffed with certified medical technologists


OSOM[®]

Trichomonas Rapid Test

Sekisui Diagnostics, LLC (AMERICAS)
 4 Hartwell Place Lexington, MA 02421
 Phone: 800 332 1042
 Fax: 800 762 6311
 Email: osomrapidtests@sekisui-dx.com
 ISO 13485 Certified

Sekisui Diagnostics (UK) Limited (INTERNATIONAL)
 Liphook Way, Allington, Maidstone, Kent ME16 0LQ, UK
 Phone: +44 1622 607800 Fax: +44 1622 607801
 Email: info@sekisui-dx.com
 ISO 9001 Certified, ISO 13485 Certified

TEST PROCEDURE / RESULTS

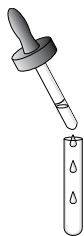


Add sample buffer to test tube.



Insert swab and mix.



Let stand 1 MINUTE.

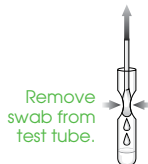

Remove swab from test tube.

Insert test stick into extracted sample.



Read results at 10 MINUTES.



POSITIVE



NEGATIVE

DESCRIPTION

OSOM[®] TRICHOMONAS RAPID TEST
 TRICHOMONAS POSITIVE CONTROL KIT

CONFIGURATION

25 TESTS

CATALOG #

181E (CE MARKED)
 182 (CE MARKED)

STORAGE

ROOM TEMP
 ROOM TEMP

