



Diagnostics Products



Table of Contents

Diagnostics Products

Goldenberg Snarecoil™* Bone Marrow Needles & Trays with Safety Components

Goldenberg Snarecoil™* Bone Marrow Needle	1
Goldenberg Snarecoil™* Bone Marrow Biopsy and Aspiration Trays with Safety Components..	1
Goldenberg Snarecoil™* Basic Bone Marrow Components	2
Basic Bone Marrow Trays with Safety Components.....	2

Monoject™ Bone Marrow Needles & Trays with Safety Components

Monoject™ I-Type Sternal-Iliac Aspiration Needle	3
Monoject™ J-Type Bone Marrow Biopsy/Aspiration Needles	3
Monoject™ Bone Marrow Biopsy Trays with Safety Components.....	4
Monoject™ Bone Marrow Aspiration Needles, Rosenthal Type	5
Monoject™ Bone Marrow Aspiration Trays, Rosenthal Type Needle, Safety Components.....	5

Argyle™ Turkel™ Safety Thoracentesis System & Trays

Argyle™ Turkel™ Safety Thoracentesis System	6
Argyle™ Turkel™ Safety Thoracentesis Procedure Trays With Safety Components	6
Argyle™ Turkel™ Safety Paracentesis Tray with Safety Components	7

Argyle™ Turkel™ Safety Pneumothorax Tray

Argyle™ Turkel™ Safety Pneumothorax Tray with Safety Components	8
---	---

Curity™ Thoracentesis Trays

Curity™ Thoracentesis Trays with Safety Components	9
--	---

Curity™ Lumbar Puncture Trays with Safety Components

Curity™ Lumbar Puncture Trays with Safety Components	10
--	----

INDEX	12
--------------------	-----------

Goldenberg Snarecoil™* Bone Marrow Needles & Trays with Safety Components



Goldenberg Snarecoil™* Bone Marrow Needle

- Features the revolutionary Snarecoil™* design that allows for the specimen to be severed and secured in the needle lumen
- Longer specimen obtained
- Less specimen fragmentation
- Stylet features a trocar-tip design for easier penetration of the cortical layer of the bone
- Snare is deployed by moving lever from “Off” to “On” position
- “In and Out” procedure - no need to “rotate” the needle to sever the specimen

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
SC5008X	8ga x 5"	10
SC4011X	11ga x 4"	10
SC6011X	11ga x 6"	10
SC2013X	13ga x 2"	10
SC3513X	13ga x 3"	10



Goldenberg Snarecoil™* Bone Marrow Biopsy and Aspiration Trays with Safety Components — Latex Free

Sterile Contents

Prep Tray

- (1) Fenestrated Drape
- (1) Towel
- (1) 3mL Applicator 2% ChloroPrep™* Wand

Procedure Tray

- (1) 6 mL Syringe with 20ga x 1 ½" Blunt Tip Filter Aspiration Needle
- (1) 20 mL Syringe and Tip Cap
- (1) 21ga x 1 ½" Magellan™ Safety Needle
- (1) 25ga x ⅝" Magellan™ Safety Needle
- (1) Probe
- (1) Probe Guide
- (1) #11 Retractable Scalpel Blade
- (2) Specimen tubes containing Edetate (EDTA) Potassium
- (10) Frosted Slides
- (1) Amber Specimen Bottle with Label
- (1) Stirring Rod
- (3) 12-Ply 3" x 3" Gauze Sponges
- (1) 1" x 3" Curity™ Adhesive Bandage
- (1) Tip Cap for Goldenberg Snarecoil™* Stylet for Luer Lock Port

With I-Type Sternal/Iliac Aspiration Needle and Bone Marrow Biopsy/Aspiration Needle

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
SC4508	16ga adjustable 3/16" to 1 7/8" Sternal/Iliac AND 8ga x 5" Bone Marrow	10
SC4511	16ga adjustable 3/16" to 1 7/8" Sternal/Iliac AND 11ga x 4" Bone Marrow	10

With Bone Marrow Biopsy/Aspiration Needle Only

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
SC3411	11ga x 4" Bone Marrow	10
SC3353	13ga x 3" Bone Marrow	10

Goldenberg Snarecoil™* Bone Marrow Needles & Trays with Safety Components

Basic Bone Marrow Trays with Safety Components — Latex Free

Sterile Contents

Prep Tray

- (1) Fenestrated Drape
- (1) Towel
- (1) 3mL Applicator 2% Chloraprep™* Wand

Procedure Tray

- (1) 6 mL Syringe with 20ga x 1 ½" Blunt Tip Filter Aspiration Needle
- (1) 20 mL Syringe
- (1) 21ga x 1 ½" Magellan™ Safety Needle
- (1) 25ga x 5/8" Magellan™ Safety Needle
- (1) Probe
- (1) Probe Guide
- (1) #11 Retractable Scalpel Blade
- (3) 12-Ply 3" x 3" Gauze Sponges
- (1) 1" x 3" Curity™ Adhesive Bandage
- (1) Tip Cap for Goldenberg Snarecoil™ Stylet for Luer Lock Port


With I-Type Sternal/Iliac Aspiration Needle and Bone Marrow Biopsy/Aspiration Needle

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
SCE4511	16ga adjustable 3/16" to 1 7/8" Sternal/Iliac AND 11ga x 4" Bone Marrow	10

With Bone Marrow Biopsy/Aspiration Needle Only

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
SCE3411	11ga x 4" Bone Marrow	10

Monoject™ Bone Marrow Needles & Trays with Safety Components

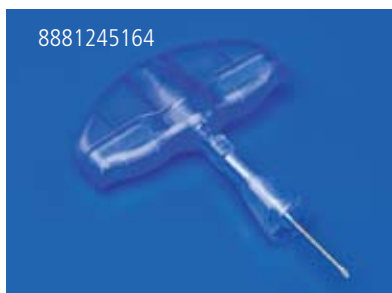

Monoject™ I-Type Sternal-Iliac Aspiration Needle

- Ergonomic handle design
- Curved finger grips
- Superior cutting bevel
- Interlocking stylet handle
- Integral luer lock fitting for luer slip or luer lock syringe
- Adjustable sternal guard for sternal applications
- Quick release sternal guard
- Sterile

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
8881245164	16ga adjustable 3/16" to 1 7/8"	10



Monoject™ J-Type Bone Marrow Biopsy/Aspiration Needles

- Ergonomic handle design
- Curved finger grips
- Superior cutting bevel
- Beveled stylet tip projects beyond needle edge
- CM markings
- Interlocking stylet handle
- Integral luer lock fitting for luer slip or luer lock syringe
- Sterile

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
8881247087	8ga x 4"	10
8881247111	11ga x 4"	10
8881247129	13ga x 2 1/2"	10
8881247137	13ga x 3 1/2"	10

Monoject™ Bone Marrow Needles & Trays with Safety Components



Monoject™ Bone Marrow Biopsy Trays with Safety Components

— **Latex Free**

Sterile Contents

Prep Tray

- (1) Fenestrated Drape
- (1) Towel
- (1) 3mL Applicator 2% Chloraprep™* Wand

Procedure Tray

- (1) 6 mL Syringe with 20ga x 1" Blunt Tip Aspiration Needle
- (1) 20 mL Syringe
- (1) 21ga x 1 ½" Magellan™ Safety Needle
- (1) 25ga x ⅝" Magellan™ Safety Needle
- (1) #11 Retractable Scalpel Blade
- (10) Frosted Slides
- (3) 12-Ply 3" x 3" Gauze Sponges
- (1) 1" x 3" Curity™ Adhesive Bandage
- (2) Specimen Tubes Containing Edetate (EDTA) Potassium
(Note: 8881847167 And 8881847993 Have 1 EDTA Tube)
- (1) Specimen Tube *without* EDTA (8881847993 Only)
- (1) Amber Specimen Bottle with Label (Not 8881847167)
- (1) Stirring Rod (Not 8881847167)
- (1) Tip Cap for Monoject™ Bone Marrow Stylet for Luer Lock Port

Tray with I-Type Sternal/Iliac Aspiration Needle and Bone Marrow Biopsy/Aspiration Needle

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
8881847019	16ga adjustable 3/16" to 1 7/8" Sternal/Iliac AND 11ga x 4" Bone Marrow	10
8881847035	16ga adjustable 3/16" to 1 7/8" Sternal/Iliac AND 8ga x 4" Bone Marrow	10

Tray with Bone Marrow Biopsy/Aspiration Needle Only

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
8881847084	8ga x 4" Bone Marrow	10
8881847118	11ga x 4" Bone Marrow	10

Tray with I-Type Sternal/Iliac Aspiration Needle Only

CODE	SHIP CASE
8881847167	16ga adjustable 3/16" to 1 7/8" Sternal/Iliac Needle 10

Bone Marrow Prep Tray Only (Does Not Include Bone Marrow Needles)

CODE	SHIP CASE
8881847993	10

Monoject™ Bone Marrow Needles & Trays with Safety Components



Monoject™ Bone Marrow Aspiration Needles, Rosenthal Type

- Plastic hub
- Traditional Bone Marrow Aspiration Needle Design
- Precision-ground and fitted needle point and stylet tip

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
8881240009	16ga x 2" Bone Marrow Needle	20
8881240017	16ga x ¾" Bone Marrow Needle	20



Monoject™ Bone Marrow Aspiration Trays with Rosenthal Type Needle with Safety Components — *Latex Free*

Sterile Contents

Prep Tray

- (1) Fenestrated Drape
- (1) Towel
- (1) 3mL Applicator 2% Chloraprep™* Wand

Procedure Tray

- (1) Bone Marrow Aspiration Needle (Rosenthal Type Metal Hub)
- (1) 6 mL Syringe with 20ga x 1 ½" Blunt Tip Filter Aspiration Needle
- (1) 20 mL Syringe
- (1) 21ga x 1 ½" Magellan™ Safety Needle
- (1) 25ga x ⅝" Magellan™ Safety Needle
- (1) Specimen Tube Containing Edetate (EDTA) Potassium
- (1) #11 Retractable Scalpel Blade
- (10) Frosted Slides
- (3) 12-Ply 3" x 3" Gauze Sponges
- (1) 1" x 3" Curity™ Adhesive Bandage

CODE	NEEDLE GAUGE AND LENGTH	SHIP CASE
8881846169	16ga x 2" Bone Marrow Needle	10

Argyle™ Turkel™ Safety Thoracentesis System & Trays



Argyle™ Turkel™ Safety Thoracentesis System

- Spring-loaded blunt safety cannula extends beyond needle tip
- One-way valve in needle housing
- Color change indicator to identify blunt cannula position
- Ball valve in catheter housing
- In-line stop cock
- Polyurethane catheter-over-needle design
- 1cm depth marks on catheter shaft
- Multiple side holes on catheter shaft
- Sterile

CODE	CATHETER SIZE	CATHETER O.D.	CATHETER LENGTH	SHIP CASE
566034	8 Fr.	2.7mm	3 ½" (9cm)	20



Argyle™ Turkel™ Safety Thoracentesis Procedure Trays With Safety Components — Latex Free

Sterile Contents

Prep Tray

- (1) CSR Wrap
- (1) Absorbent Towel
- (1) 3mL Applicator 2% Chloraprep™* Wand
- (1) Fenestrated Drape
- (6) 12-Ply 3" x 3" Gauze Sponges

Procedure Tray

- (1) 8Fr. Turkel™ Safety Thoracentesis System
- (1) 6 mL, 12 mL, 60 mL Monoject™ Syringe
- (1) 25ga x 5/8" Magellan™ Safety Needle
- (1) 21ga x 1 ½" Magellan™ Safety Needle
- (1) 20ga x 1 ½" Blunt Tip Filter Aspiration Needle
- (1) #11 Retractable Scalpel Blade
- (3) 8mL Specimen Vials
- (1) Tip Cap
- (1) 1" x 3" Curity™ Adhesive Bandage
- (1) Drainage Collection Pouch:
 - Automatic aspiration/purge valve and 36" tubing set with set-up diagram
 - Standard 36" tubing set and roller clamp with set-up diagram
 - 2000mL drainage collection bag

CODE	CATHETER SIZE	CATHETER O.D.	CATHETER LENGTH	SHIP CASE
566059	8 Fr.	2.7mm	3 ½" (9cm)	5
566075	8 Fr.	2.7mm	5" (12.5cm)	5

Argyle™ Turkel™ Safety Paracentesis Tray



Argyle™ Turkel™ Safety Paracentesis Tray with Safety Components

— Latex Free

Turkel™ Paracentesis Safety Needle

- Blunt, multi-side holed, spring-loaded inner cannula extends beyond needle bevel; protects underlying organs from inadvertent penetration
- Color change indicator, located in the needle hub provides direct visualization of the position of the blunt cannula, as the needle moves through the abdominal wall and into the peritoneal space

Paracentesis Catheter

- Polyurethane, radiopaque catheter material provides strength, softness flexibility, and resistance to thrombosis
- Multiple side holes in distal catheter shaft increase surface area for drainage
- Depth marks on catheter shaft provide assessment of catheter location during procedure
- Tuohy-Borst adapter provides a way to reinsert the safety needle for multiple insertion capability

Paracentesis Procedure Tray

- Preassembled disposable tray saves procedural set-up time; assures availability of a complete, sterile tray
- Two-tiered tray design separates prep and procedure components; easy procedure planning
- Sequential arrangement of components reduces procedure time

Sterile Contents

Prep Tray

- (1) CSR Wrap
- (1) Absorbent Towel
- (1) 3mL Applicator 2% Chloraprep™* Wand
- (1) Fenestrated Drape
- (6) 12-Ply 3" x 3" Gauze Sponges

Procedure Tray

- (1) 8 Fr. Turkel™ Safety Paracentesis System
- (1) 6 mL, 12 mL, 60 mL Monoject™ Syringe
- (1) 25ga x 5/8" Magellan™ Safety Needle
- (1) 21ga x 1 1/2" Magellan™ Safety Needle
- (1) 20ga x 1 1/2" Blunt Tip Filter Aspiration Needle
- (1) 16ga x 1 1/2" Plastic Hub Needle
- (1) #11 Retractable Scalpel Blade
- (3) 8 mL Specimen Vials
- (1) 13 x 75mm, 5mL EDTA Blood Collection Tube
- (1) Tip Cap
- (1) 1" x 3" Curity™ Adhesive Bandage
- (1) Drainage Collection Pouch:
 - (1) Vented 2000ml Drainage Bag
 - (1) 36" Tubing Set with Roller Clamp
 - (1) 36" Tubing Set with Automatic Aspiration/Purge Valve

CODE	CATHETER SIZE	CATHETER O.D.	CATHETER LENGTH	SHIP CASE
8888568006	8 Fr.	2.7mm	4 3/4" (12cm)	5

Argyle™ Turkel™ Safety Pneumothorax Tray



Argyle™ Turkel™ Safety Pneumothorax Tray with Safety Components

Safety Pneumothorax Needle

- Blunt, multi-side holed, spring-loaded inner cannula extends beyond needle bevel; protects underlying tissue (lung) from further, inadvertent penetration
- One-way duck bill valve in needle hub housing allows air and fluid to escape from the needle yet prevents atmospheric air from entering into the pleural space to prevent accidental pneumothorax during procedure
- Color change indicator in needle hub housing provides direct visualization of the position of the cannula as the needle moves through the chest wall and into the pleural space

Pneumothorax Catheter

- Polyurethane catheter material provides strength, softness, flexibility and resistance to thrombosis
- Multiple side holes in distal catheter shaft increases surface area for drainage
- Depth marks on catheter shaft provide assessment of catheter location during procedure

Pneumothorax Procedure Tray

- Preassembled disposable tray saves procedural setup time; assures availability of a complete, sterile tray
- Two-tiered procedure tray design separates prep and procedure components; easy procedure planning
- Sequential arrangement of components reduces procedure time

Sterile Contents

Prep Tray

- (1) CSR Wrap
- (1) Absorbent Towel
- (1) 3mL Applicator 2% Chloraprep™* Wand
- (1) Fenestrated Drape
- (6) 12-Ply 3" x 3" Gauze Sponges

Procedure Tray

- (1) 8 Fr. Turkel™ Safety Pneumothorax System (Percutaneous Thoracic Catheter)
- (1) 6 mL, 12 mL Monoject™ Syringe
- (1) 25ga x 5/8" Magellan™ Safety Needle
- (1) 21ga x 1 1/2" Magellan™ Safety Needle
- (1) 20ga x 1 1/2" Blunt Tip Filter Aspiration Needle
- (1) #11 Retractable Scalpel Blade
- (1) Heimlich Valve*
- (1) 36" Tubing Set with Funnel Connector
- (1) Tube Adapter*
- (1) One-Way Stopcock
- (1) 1" x 3" Curity™ Adhesive Bandage
- (1) Silk Suture with Cutting Needle

* All items in tray are LATEX-FREE except for the Heimlich Valve and Tube Adapter, which contain latex.

CODE	CATHETER SIZE	CATHETER O.D.	CATHETER LENGTH	SHIP CASE
8888567032	8 Fr.	2.7mm	3 1/2" (9cm)	5

Curity™ Thoracentesis Trays



Curity™ Thoracentesis Trays with Safety Components — *Latex Free*

- Aspirating needle is equipped with with a unique depth clamp that adjusts easily to the desired level of insertion and then locks into place, reducing the risk of over-penetration of the needle into the chest cavity
- Once assembled, the closed system allows pleural fluid to be aspirated, and collected with reduced risk of external contamination

Sterile Contents

Prep Tray

- (1) Underpad
- (1) 3mL Applicator 2% ChloroPrep™* Wand
- (1) Removable Prep Tray

Procedure Tray

- (3) 12-Ply 3" x 3" Gauze Sponges
- (1) 25ga x 5/8" Magellan™ Safety Needle
- (1) 21ga x 1 1/2" Magellan™ Safety Needle
- (1) Plastic Syringe, 5mL Luer Slip
- (1) Plastic Syringe, 50mL Luer Lock
- (1) Stopcock/Side Arm Assembly
- (3) 8 mL Specimen Vials & Caps
- (1) Collection Bag, 1500mL with Quick Disconnect
- (1) Fenestrated Drape
- (1) 1" x 3" Curity™ Adhesive Bandage

CODE	DESCRIPTION	SHIP CASE
5014	Includes 16ga x 3 1/2" Aspirating Needle with Adjustment Clamp	10
5016	Includes 17ga x 6" Aspirating Needle with 14ga Teflon and Adjustable Depth Guard	10

Curity™ Amniocentesis Tray



Curity™ Amniocentesis Tray — *Latex Free*

- Exclusive Amniocentesis needle for safe and effective aspiration
- Disposable trays contain all the components necessary for aspiration of amniotic fluid

Sterile Contents

- (1) 20ga x 3 ½" Needle with Stylet
- (1) CSR Wrap
- (1) Component Tray with Metric Measures
- (2) Sponge Sticks
- (3) Povidone Iodine Swabsticks
- (4) 12-Ply 3" x 3" Gauze Sponges
- (1) Fenestrated Drape
- (1) 22ga x 1 ½" Magellan™ Safety Needle
- (1) 25ga x 1" Magellan™ Safety Needle
- (1) 5mL Luer Lock Syringe
- (1) 20mL Luer Lock Syringe
- (1) Amber Pouch
- (3) Sample Collection Tubes with Caps
- (2) Sample Labels
- (1) 1" x 3" Curity™ Adhesive Bandage

CODE	DESCRIPTION	SHIP CASE
97051	20ga x 3 ½" Amniocentesis Needle	10

Curity™ Lumbar Puncture Trays with Safety Components



Curity™ Lumbar Puncture Trays with Safety Components — *Latex Free*

- Easy attachment of the manometer and sample vials to the adapter facilitates quick assembly and eliminates spillage of cerebral spinal fluid
- Exclusive automatic stopcock contains self-sealing valves that open and close automatically to ensure an aseptic process
- Extension line reduces manometer torque on lumbar puncture needle, improving patient comfort
- Easy one-handed manipulation of sample vials

Sterile Contents

Prep Tray

- (1) Underpad
- (3) 12-Ply 3" x 3" Gauze Sponges
- (1) Removable Prep Tray
- (2) Prep Sponge Applicator Assembly

Procedure Tray

- (1) 25ga x 5/8" Magellan™ Safety Needle
- (1) 21ga x 1 1/2" Magellan™ Safety Needle (Adult + Pediatric Only)
- (1) 3mL Plastic Syringe, Luer Slip
- (1) Automatic Adapter
- (1) Two-Piece Manometer, 56cm (Adult + Pediatric Only)
- (4) Sample Vials, with Caps
- (1) Fenestrated Drape
- (1) 1" x 3" Curity™ Adhesive Bandage
- (1) Extension Tubing

CODE	DESCRIPTION	WITH SPINAL NEEDLE	SHIP CASE
<i>Adult Lumbar Puncture Trays, Sterile</i>			
1017	Adult Lumbar Puncture Tray	18ga x 3 1/2", with Stylet	20
1032	Adult Lumbar Puncture Tray	20ga x 3 1/2", with Stylet	20
<i>Pediatric Lumbar Puncture Trays, Sterile</i>			
1048	Pediatric Lumbar Puncture Tray	22ga x 2 1/2", with Stylet	10
<i>Infant Lumbar Puncture Trays, Sterile</i>			
7039	Infant Lumbar Puncture Tray	22ga x 1 1/2", with Stylet	10
7500	Infant Lumbar Puncture Tray, with Manometer 29cm	22ga x 1 1/2", with Stylet	10

INDEX

CODE	PRODUCT	PAGE
1017	Curity Lumbar Puncture Safety Tray	11
1032	Curity Lumbar Puncture Safety Tray	11
1048	Curity Lumbar Puncture Safety Tray	11
5014	Curity Thoracentesis Safety Tray	9
5016	Curity Thoracentesis Safety Tray	9
566034	Argyle Turkel Safety Thoracentesis System	6
566059	Argyle Turkel Safety Thoracentesis Tray	6
566075	Argyle Turkel Safety Thoracentesis Tray	6
7039	Curity Lumbar Puncture Safety Tray	11
7500	Curity Lumbar Puncture Safety Tray	11
8881240009	Monoject Bone Marrow Aspiration Needle	5
8881240017	Monoject Bone Marrow Aspiration Needle	5
8881245164	Monoject I-Type Sternal-Iliac Aspiration Needle	3
8881247087	Monoject J-Type Bone Marrow Biopsy/Aspiration Needles	3
8881247111	Monoject J-Type Bone Marrow Biopsy/Aspiration Needles	3
8881247129	Monoject J-Type Bone Marrow Biopsy/Aspiration Needles	3
8881247137	Monoject J-Type Bone Marrow Biopsy/Aspiration Needles	3
8881846169	Monoject Bone Marrow Aspiration Safety Tray	5
8881847019	Monoject Bone Marrow Biopsy Safety Trays	4
8881847035	Monoject Bone Marrow Biopsy Safety Trays	4
8881847084	Monoject Bone Marrow Biopsy Safety Trays	4
8881847118	Monoject Bone Marrow Biopsy Safety Trays	4
8881847167	Monoject Bone Marrow Biopsy Safety Trays	4
8881847993	Monoject Bone Marrow Biopsy Safety Trays	4
8888567032	Argyle Turkel Safety Pneumothorax Tray	8
8888568006	Argyle Turkel Safety Paracentesis Tray	7
97051	Curity Amniocentesis Tray	10
SC2013X	Goldenberg Snarecoil Bone Marrow Needle	1
SC3353	Goldenberg Snarecoil Safety Tray	1
SC3411	Goldenberg Snarecoil Safety Tray	1
SC3513X	Goldenberg Snarecoil Bone Marrow Needle	1
SC4011X	Goldenberg Snarecoil Bone Marrow Needle	1
SC4508	Goldenberg Snarecoil Safety Tray	1
SC4511	Goldenberg Snarecoil Safety Tray	1
SC5008X	Goldenberg Snarecoil Bone Marrow Needle	1
SC6011X	Goldenberg Snarecoil Bone Marrow Needle	1
SCE3411	Goldenberg Snarecoil Basic Safety Tray	2
SCE4511	Goldenberg Snarecoil Basic Safety Tray	2

The products shown in this catalog may not be licensed or registered in all parts of the world. Please contact your local sales representative for availability.

TM* CHLORAPREP is a trademark of Allegiance Corp.

TM* GOLDENBERG SNARECOIL is a trademark of Ranfac Corp.

TM*TEFLON is a trademark of E.I. DuPont de Nemours & Co.

COVIDIEN, COVIDIEN with Logo and TM marked brands are trademarks of Covidien AG or an affiliate.

©2009 Covidien. All rights reserved.



H6200DIAG1M-0809