

Hot/Cold Therapy & Clinic Supplies





# Introduction

Chattanooga is the world's largest manufacturer of rehabilitation equipment for treating musculoskeletal, neurological and soft tissue disorders. For over six decades, Chattanooga has set a benchmark for leadership, reliability and excellence. Our products contribute to better treatment outcomes in hospitals, clinics and home settings worldwide. We lead the rehabilitation industry by example, and through continuous innovation, provide real world solutions for clinicians and their patients.

As the world's leading brand in rehabilitation, Chattanooga is sold in over 80 countries worldwide. Chattanooga distributors and customers are supported by a highly trained and seasoned team of personnel who oversee new product development, product training and education, order placement, product delivery, and after sales service in addition to providing dedicated sales and marketing support.

Chattanooga is a brand of DJO Global Inc. DJO Global provides solutions for musculoskeletal health, vascular health and pain management. Our products help prevent injuries or rehabilitate after surgery, injury or from degenerative disease, enabling patients to regain or maintain their natural motion.

Visit www.DJOglobal.com/chattanooga

















# Contents

Pedal Exerciser

Stabilizer

Saunders Total Back

Chattanooga Evaluation Devices

Hydrocollator Mobile Heating Units	2	Baseline Evaluation Devices	1
Hydrocollator Stationary Heating Units	3	Sully Shoulder Support	1
Hydrocollator HotPacs	4	Gels & Lotions	1
Hydrocollator Terry Covers	5	Electrodes & Accessories	1
ColPac Chilling Unit	6	Education Manuals	Ź
ColPac, Boo-Boo Pac, Flexipac, Sport-Pac	7		
Theratherm Digital Moist Heat	8		
Carts & Stools	9		

10

10

11

12



Hot/Cold Therapy and Clinic Supplies

# **Ordering Information** Item# Description M-2 Mobile Hydrocollator M2 Mobile - Includes 1 2402 2 Standard (1006) HotPacs M2 Mobile - Includes 3 Standard (1006), 402-2 3 Oversize (1004), and 3 Cervical (1002) HotPacs M2 Mobile - Includes 4 Standard (1006), 2402-3 2 Oversize (1004), and 4 Cervical (1002) HotPacs M2 Mobile - Includes 4 Standard (1006), and 2402-4 4 Oversize (1004) M-4 Mobile Hydrocollator M4 Mobile - Includes 24 Standard (1006) HotPacs M4 Mobile - Includes 8 Oversize (1004), and 8 Cervical (1002) HotPacs M4 Mobile - Includes 6 Standard (1006), 6 Oversize (1004), and 6 Cervical (1002) HotPacs SS-2 Mobile Hydrocollator SS2 Mobile - Includes 8 Standard (1006) HotPacs SS2 Mobile - Includes 4 Standard (1006), and

# SS2 Mobile - Includes 2 Standard (1006), 2 Oversize (1004), and 2 Cervical (1002) HotPacs **Technical Specifications** M-2 Mobile Hydrocollator Mains Power: 110~120 V, 50/60 Hz Power Consumption: 1000 W 48 lb (22 kg) /Shipping - 93 lb (42 kg) Weight: Dimensions\*: 27" x 16" x 33" (69 x 41 x 84 cm) Electrical Safety Class: Class 1, Type B Conforms to UL 60601-1, certified to Safety Tests: Can/CSA C222 No. 601.1 M-4 Mobile Hydrocollator

4 Cervical (1002), HotPacs

Mains Power:	110~120 V, 50/60 Hz
Power Consumption:	1500 W
Weight:	132 lb (60 kg) /Shipping - 160 lb (72 kg)
Dimensions*:	35" x 20" x 33" (89 x 51 x 84 cm)
Electrical Safety Class:	Class 1, Type B
Safety Tests:	Conforms to UL 60601-1, certified to Can/CSA C222 No. 601.1

Electrical surcey class.	ciass i, Type b
Safety Tests:	Conforms to UL 60601-1, certified to Can/CSA C222 No. 601.1
SS-2 Mobile Hydrocollato	or
Mains Power:	110~120 V, 50/60 Hz
Power Consumption:	1000 W
Weight:	66 lb (30 kg) /Shipping - 80 lb (36 kg)
Dimensions*:	21" x 16" x 33" (53 x 41 x 84 cm)
Electrical Safety Class:	Class 1, Type B
Safety Tests:	Conforms to UL 60601-1, certified to Can/CSA C222 No. 601.1
* Dimensions are listed a	ıs L x W x H

# Hydrocollator® Mobile Heating Units



Our benchmark product is the standard all others can only hope to match. Durable and easy to maintain, these high-quality stainless steel units give you a constant supply of temperature-consistent HotPacs.

### **Features**

#### **Hydrocollator Mobile Heating Units**

- High-quality stainless steel
- · Thermostatically controlled temperature
- 3" (8 cm) swivel-type rubber casters for silent, friction-free movement
- Easy maintenance. Simple to fill and drain. No plumbing required
- One-year warranty

\* Approximate with HotPacs

				- ''	prominate mairrioti des
	M-2	M-4	SS-2	E-2	E-1/E-1 MHT
Tank Capacity	14 gal (52 L)	27 gal (102 L)	10 gal (37 L)	8 gal (37 L)	3 gal (37 L)
Temperature Range	160° - 165°F (71° - 74°C)				
Thermal Cut-Out Temp.	180° - 185°F (82° - 85°C)				
Temperature Accuracy	+/-10%	+/- 10%	+/- 10%	+/-10%	+/- 10%
Heat up Time to 160° F (to 70°C)	6 Hours	8 Hours	4 Hours	4 Hours	2 Hours
Cool Down Time from 160° F (from 70°C)	3 Hours	4 Hours	4 Hours	2 Hours	1 Hour
Fiberglass Insulation	Yes	Yes	Yes	No	No

# Hydrocollator® Stationary Heating Units



# Hydrocollator Accessories



### **Wall Mounted Towel Rack**

Six hooks extend 13" (33 cm) from wall. Measures 2" x 33" (5 cm x 84 cm).



#### Side Table Rack

Stainless steel work surface and two wire shelves add efficiency to units. Use with M-4, M-2 and SS-2.



### Mini Non-Electric Lotion Warmer

Hangs on the front of the Hydrocollator E-1 or E-2 heating unit and uses the unit's warmth to keep lotion or gel at a comfortable temperature. Bottles are not included.



### **TM-1 Electric Lotion Warmer**

Easily attaches to any treatment table or wall to keep massage lotion or ultrasound gel at a comfortable temperature and within reach. Bottles are not included.

Ordering	Information	
Item#	Description	
E-1 Stationa	ary Hydrocollator	
2102	E1 Stationary - Includes 4 Standard (1006) HotPacs	
2102-2	E1 Stationary - Includes 3 Standard (1006) HotPacs and 1 Cervical (1002) HotPac	
2102-2	E1 Stationary - Includes 2 Standard (1006) HotPacs and 2 Cervical (1002) HotPacs	
E-2 Station	ary Hydrocollator	
2802	E1 Stationary - Includes 2 Oversize (1004), 3 Standard (1006) HotPacs and 1 Cervical (1002) HotPac	
Hydrocollat	tor Accessories	
4016	Wall Mounted Towel Rack	
4010	Side Table Rack	

Extra Shelf for Side Table Rack (not shown)

Mini Non-Electric Lotion Warmer

Dial Thermometer (not shown)

TM-1 Electric Warmer

# **Technical Specifications**

E 2 Stationary Hydrocollator

4230

4012

4228

E-1 Stationary Hydrocol	lator
Mains Power:	110~120 V, 50/60 Hz
Power Consumption:	1000 W
Weight:	15 lb (7 kg) /Shipping - 21 lb (10 kg)
Dimensions*:	13" x 8" x 16" (33 x 20 x 41 cm)
Electrical Safety Class:	Class 1, Type B
Safety Tests:	Conforms to UL 60601-1, certified to Can/CSA C222 No. 601.1

Mains Power: 110~120 V, 50/60 Hz  Power Consumption: 1000 W  Weight: 29 lb (13 kg) /Shipping - 46 lb (21 kg)  Dimensions*: 15" x 13" x 20" (38 x 33 x 51 cm)  Electrical Safety Class: Class 1, Type B  Safety Tests: Conforms to UL 60601-1, certified to Can/CSA C222 No. 601.1	E-2 Stationary Hydrocollator			
Weight: 29 lb (13 kg) /Shipping - 46 lb (21 kg)  Dimensions*: 15" x 13" x 20" (38 x 33 x 51 cm)  Electrical Safety Class: Class 1, Type B  Safety Tests: Conforms to UL 60601-1, certified to	Mains Power:	110~120 V, 50/60 Hz		
Dimensions*: 15" x 13" x 20" (38 x 33 x 51 cm)  Electrical Safety Class: Class 1, Type B  Safety Tests: Conforms to UL 60601-1, certified to	Power Consumption:	1000 W		
Electrical Safety Class: Class 1, Type B  Safety Tests: Conforms to UL 60601-1, certified to	Weight:	29 lb (13 kg) /Shipping - 46 lb (21 kg)		
Conforms to UL 60601-1, certified to	Dimensions*:	15" x 13" x 20" (38 x 33 x 51 cm)		
CATALY LACTS.	Electrical Safety Class:	Class 1, Type B		
	Safety Tests:	•		

\* Dimensions are listed as L x W x H

Ordering Information		
Item#	Description	
Hydrocolla	ator HotPac - Moist Heat Hot Pack	
A 1008	10" x 18" (25 cm x 46 cm)	
B 1010	10" x 24" (25 cm x 61 cm)	
C 1004	Oversize 15" x 24" (38 cm x 61 cm)	
D 1002	Neck Contour 24" long (61 cm)	
E 1014	Half Size 5" x 12" (13 cm x 30 cm)	
F 1006	Standard 10" x 12" (25 cm x 30 cm)	
G 1012	Knee-Shoulder 10" x 20" (25 cm x 50 cm)	

# Hydrocollator® Hotpacs™





The first choice for physical therapists, physicians, athletic trainers, and chiropractors for over 60 years.

### Features

# The Original Hydrocollator® Moist Heat HotPac

Chattanooga is the inventor of the original Hydrocollator HotPac and proud to still be producing this all natural product today. The benefits of moist heat therapy over dry heat are well established for pain relief due to arthritis, myofibrosis, back pain, muscle strains and spasms.

#### Hydrocollator HotPacs

- Seven HotPac sizes designed to fit every body contour
- Provides up to 30 minutes of deep, soothing, moist heat
- Brown tabs on most HotPacs for easy rotation in busy clinics
- One HotPac provides hundreds of cost effective treatments
- One-year warranty







# Hydrocollator® Terry Covers





### **Features**

# **All-Terry Covers**

- For high-volume institutional use or therapy at home
- Hook and loop closures to prevent unwrapping or slipping
- · Washable and reusable
- · One-year warranty

# **Foam-Filled Terry Covers**

- For home or light clinic use
- Polyurethane foam layer between two layers of terry
- Available with or without pockets
- Washable and reusable
- One-year warranty





Description
vers
Standard 19" x 27" (48 cm x 69 cm)
Oversize 24" x 36" (61 cm x 91 cm)
Neck Contour 25" x 18" (64 cm x 46 cm)
Terry Covers
Standard 19" x 27" (48 cm x 69 cm)
Oversize 24" x 36" (61 cm x 91 cm)
Neck Contour 25" x 16" (64 cm x 41 cm)
ket Terry Covers
Standard 17" x 27" (43 cm x 69 cm)
Oversize 24" x 36" (61 cm x 91 cm)
or HotPac Set in Retail Packaging
Standard Set in Retail Package includes one Standard HotPac and one Foam Filled Terry Cover

Neck Set in Retail Package includes one Neck HotPac and one Foam

Filled Terry Cover

1062

# Ordering Information

Item# Description

# ColPac Chilling Unit Model C-5

3102 C-5 ColPac Chilling Unit Includes 6 Standard (1500) and 6 Half Size (1506) ColPacs

#### Standard Accessories

Side Table Rack See page 4 for Description and image

4230 Extra Shelf for Side Table Rack (not shown)

#### ColPac Freezer

ColPac Freezer 21 ¼" x 26" x 33½" 90910 (53 cm x 66 cm x 86 cm) Includes 12 Standard ColPacs (1500) Shipping weight 150 lbs (68 kg)

### **Technical Specifications**

ColPac Chilling Unit Model C-5			
Mains Power:	115 VAC~ 60 Hz		
Current:	4A		
Weight:	165 lb (74 kg) Shipping Weight - 190 lb (86 kg)		
Dimensions*:	27" x 16" x 33" (69 x 41 x 84 cm)		
Electrical Safety Class:	Туре В		
Safety Tests:	UL 544 (U.S. Only)		

\* Dimensions are listed as LxWxH

# ColPac® Chilling Unit





Large capacity and fast chilling time means your patients will never have to wait for a ColPac. The ColPac Chilling Units help you keep pace throughout a busy day. No plumbing required, simply plug into a standard electrical outlet.

#### **Features**

#### **ColPac Chilling Unit**

- Closed-cell foam insulation for efficiency
- Stainless steel, welded construction for long life
- Heavy-duty compressor for more efficient cooling
- Specially designed interior chilling coils cool packs faster
- 3-inch (8 cm) swivel-type rubber casters for silent, friction-free movement of mobile units
- Drain valve for easy cleaning and defrosting
- · One-year warranty

# ColPac Freezer

- Five cubic feet of spacious capacity
- Four durable stainless steel shelves
- · Comes complete with trays and cup holders
- Shipped with 12 standard size ColPacs
- UL listed and cUL listed
- · One-year warranty





# ColPac®, Boo-Boo Pac™, Flexipac®, Sport-Pac™

Seven different sizes of blue vinyl and four sizes of black polyurethane ensure effective cold therapy for any body area or shape.



#### **Features**

### The Original ColPac

- · Designed for heavy clinical use.
- · Effective pain relief
- Filled with a non-toxic silica gel that remains pliable throughout treatment
- Delivers up to 30 minutes of soothing relief from acute pain, swelling and fever
- Latex free
- One-year warranty

# boo-boo pac

- Fun, bear-shaped design especially made for kids
- Ideal for pediatric applications
- Filled with non-toxic silica gel
- Delivers up to 30 minutes of soothing relief from acute pain, swelling and fever
- Shipped in a display box containing 12 packs

#### FlexiPac

- · Freezable to provide penetrating cold therapy
- Microwavable for soothing heat
- Flexible gel conforms to any body contour at any temperature

#### Sport-Pac

- · Soccer ball design is ideal for youth sport applications
- Easily transferred from the freezer to the team cooler for sideline use
- Delivers up to 30 minutes of effective cold therapy for acute pain due to strains, sprains and bruises
- Shipped in a display box containing 10 packs

Ordering Information		
Item#	Description	
Blue Vinyl C	olPac	
A 1512	Oversize 11" x 21" (28 cm x 53 cm)	
B 1500	Standard Size 11" x 14" (28 cm x 36 cm)	
C 1506	Half Size 7½" x 11" (19 cm x 28 cm)	
D 1508	Neck Contour 23" long (58 cm)	
E 1502	3" x 11" (8 cm x 28 cm)	
F 1504	Quarter Size 5½" x 7½" (14 cm x 19 cm)	
G 1510	Eye (Specially shaped to cover closed eyelids)	
Black Polyu	rethane ColPac	
H 1554	Neck Contour 21" long (53 cm)	
I 1556	Oversize 12½" x 18½" (32 cm x 47 cm)	
J 1552	Standard Size 10" x 13½" (25 cm x 34 cm)	
K 1562	Half Size 6½" x 11" (17 cm x 28 cm)	
L boo-boo pac		
1534	Royal Blue boo-boo pac (12 per case) Shipped in a display box containing 12 packs. NOTE: Must be purchased in quantities of 12 each	
M FlexiPac	Hot & Cold Compresses	
4026	5" x 6" (13 cm x 15 cm) 48/case	

# 4020 5" x 10" (13 cm x 25 cm) 24/case 4029 8" x 14" (20 cm x 35 cm) 12/case 5" x 10" (13 cm x 25 cm) 25/case

Chiropractically Yours (not shown)

# N Sport-Pac

4032

Sport-Pac (10 per case)
Shipped in a display box containing 10 packs.
NOTE: Must be purchased in quantities of 10 each

Ordering Information	
Item#	Description
Therather	m Digital Moist Heat Pack
1030	Small 7" x 15" (18 cm x 38 cm)
1031	Medium 14" x 14" (36 cm x 36 cm)
1032	Standard 14" x 27" (36 cm x 68 cm)
1033	Shoulder/Neck 23" x 20" (58 cm x 50 cm)
Replaceme	ent Flannel Covers
13330	Small 7" x 15" (18 cm x 38 cm)
13331	Medium 14" x 14" (36 cm x 36 cm)
13332	Standard 14" x 27" (36 cm x 68 cm)
13343	Shoulder/Neck 23" x 20" (58 cm x 50 cm)

# Theratherm® Digital Moist Heat Pack



Theratherm Digital Moist Heat Pack provides intense, moist heat therapy for temporary relief from pain associated with arthritis, back pain caused by muscle spasm and inflammation caused by strain and tension. Its flannel cover draws moisture from the air and as the pack's temperature rises, the moisture is released onto the skin, speeding relief to your patients. Its unique, easy-to-use digital control is convenient and safe.

# **Features**

# Theratherm Digital Moist Heat Pack

- Convenient preset time and temperature
- Digital hand control for operator to program and monitor temperature and/or treatment time
- Lockout mode protects against inconsistent setting changes throughout treatment sessions
- Moist heat therapy temperatures range from 88° to 166° F (32° to 74° C)
- Treatment sessions up to 60 minutes
- Specially designed flannel cover for efficient moisture retention
- · UL and cUL classified
- One-year warranty









Medium Standard Shoulder/Neck

# Stainless Steel Carts





### **Features**

#### **Stainless Steel Carts**

Versatile, functional and mobile, these sturdy carts go wherever you need them.

# Pneumatic Therapy Stools



### **Features**

# **Pneumatic Therapy Stools**

A high-quality, low-priced therapy stool with round, 3" (8 cm) thick foam upholstered swivel seat. Available with or without back in the four colors listed below. A convenient gas spring makes raising and lowering the stool simple 18" - 22" (46 cm - 56 cm). Maximum weight 350 lbs (158 kg).

Ordering Information		
Item#	Description	
Stainless Sta	eel Carts	
A 4002	Model MB – Specially designed to hold the SS, E-2 and E-1 heating units 20" H x 19" W x 17" D (51 cm x 48 cm x 43 cm)	
B 4004	Model MB-T – Accommodates stationary heating, ultrasound and electrotherapy units. 30" H x 19" W x 17" D (76 cm x 48 cm x 43 cm)	
C 4018	Model MB-TD – Accommodates stationary heating and electrotherapy units. Drawer provides convenient storage for clinic supplies. 30" H x 19" W x 17" D (76 cm x 48 cm x 43 cm)	
Pneumatic Stools with back		

//064	Imperial Blue
77065	Burgundy
77066	Dove Gray
77067	Black
Pneumatic S	tools without back

Pneumatic Stools without back		
77060	Imperial Blue	
77061	Burgundy	
77062	Dove Gray	
77063	Black	

Ordering Information		
Item#	Description	
Chattanooga S	tandard Exerciser™	
18000	Retail Package	
18020	Standard Package	
Saunders® To	otal Back Exercise Device	
640243-002	Saunders Total Back Black	
640243-003	Saunders Total Back Imperial Blue	
640243-005	Saunders Total Back Burgundy	

# Chattanooga Standard Exerciser™



#### **Features**

#### Chattanooga Standard Exerciser™

- Perfect for upper or lower body exercises
- Extremely lightweight: 5 lbs (2.25 kg)
- Wrap-around foot straps
- Comes completely assembled
- Retail packaging available
- Shipping dimensions: 9" x 16" x 45"
- One-year warranty



# Saunders® Total Back Exercise Device



### **Features**

# Saunders® Total Back Exercise Device

Promotes strength, flexibility and coordination through controlled functional exercise

- Pelvic stabilization allows isolation of lumbar extensors in any range of motion
- · Adjusting the table angle allows either mild exercise in a pain free range, or intense conditioning
- Angle of recline reduces the load on the spine by 20% to 30% and allows for spinal positions in flexion, extension or side lying though a range of 110°
- "Side sit-ups" are great for the obliques, transversus abdominus, multifundus, quadratus lumborum and erector spinae muscles
- Includes a full-size wall poster feature 28 numbered exercises

# Stabilizer™





### **Features**

### Stabilizer™

The Stabilizer Pressure Biofeedback unit gives valuable information to ensure quality and precision in exercise performance and muscle testing. It allows the clinician and patient to determine if the patient is able to selectively isolate and maintain contractions of the cervical or lumbopelvic core stabilization muscles. The measuring range is 0-200 mmHg analog pressure with an accuracy of +/- 3 mnHz pressure.

- Use with exercise to protect and stabilize joints
- Prevent and treat low back and neck pain
- Improve core stabilization exercises
- Enhance home exercise programs



# Ordering Information Item# Description Stabilizer™ Stabilizer Unit 9296 Comes complete with display gauge, pressure cell and instruction booklet.

# Ordering Information

Description Item#

### Chattanooga Hydraulic Hand Dynamometer

43050 200 lbs (91 kg) Dial Gauge

#### Hydraulic Pinch Gauge

43051 50 lbs (23 kg) Dial Gauge

### Skinfold Caliper

43063 Skinfold Caliper

# Chattanooga Hydraulic Hand Dynamometer

#### **Features**

#### Hydraulic Hand Dynamometer

- · Gives accurate grip strength readings without being able to "feel" the handle move
- · Dynamometer accurately measures grip strength up to 200 lbs (91 kg)
- · Adjustable handle with five positions to accommodate any hand size
- · Indicator remains at maximum reading until reset
- · Includes carrying case
- One-year warranty



# Chattanooga Hydraulic Pinch Gauge

#### **Features**

### **Hydraulic Pinch Gauge**

- Reliably and accurately register up to 50 lbs (23 kg)
- Measures all pinch tests (tip, key and palmer)
- · Includes carrying case
- One-year warranty



# Chattanooga Skinfold Caliper

### **Features**

### **Skinfold Caliper**

- · Constructed of lightweight aluminum
- · Large 60 mm scale is easy to read
- · Measurement arms and "floating tips" assure accurate and reproducible results
- · Includes carrying case, booklet and body fat percentage tables
- One-year warranty



# Chattanooga 3-Piece Hand Evaluation Set

#### **Features**

#### 3-Piece Hand Evaluation Set

- · All the instruments necessary to evaluate strength and range of motion for hand and finger
- · Hydraulic Hand Dynamometer
- Hydraulic Pinch Gauge (50 lbs/23 kg)
- · Stainless steel finger goniometer
- · Includes carrying case
- One-year warranty



### Ordering Information Description Item# 3-Piece Hand Evaluation Set 43055 200 lbs (91 kg) Dial Gauge Dynamometer 7-Piece Hand Evaluation Set 200 lbs (91 kg) Dial Gauge Dynamometer 30 lbs (14 kg) Pinch Gauge Goniometers A 43058 Plastic 360° ISOM 12" (30 cm) B 43059 Plastic 360° ISOM 8" (20 cm) C 43060 Plastic 360° ISOM 6" (15 cm) Stainless Steel Finger and Small Joint 6" (15 cm) D 43057

# Chattanooga 7-Piece Hand Evaluation Set

#### **Features**

#### 7-Piece Hand Evaluation Set

- · All instruments needed for a full hand evaluation
- Hydraulic Hand Dynamometer
- · 2-point discriminator with 3rd point
- · Finger circumference gauge and finger motion gauge
- Wartenburg pinwheel
- · Hydraulic pinch gauge registering up to 50 lbs (23 kg)
- · Stainless steel goniometer
- Digit-O-Meter
- · Includes carrying case
- One-year warranty




# Chattanooga Goniometers

### **Features**

### **Goniometers**

- The plastic 360° scale ISOM (International Standards of Measurement) goniometers (A, B, C) have a head with three scales calibrated for use with the ISOM system. The clear plastic permits observation of the joint's axis of motion and range of motion
- The stainless steel finger and small joint goniometer (E) measures the range of motion of finger joints (metacarpophalangeal and interphalangeal) and other small joints including toes. The head has two opposing 180° scales, which are marked in 5° increments
- · All goniometers have a linear scale in inches and centimeters
- One-year warranty

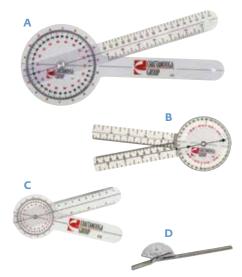

# Ordering Information

Item# Description

#### Digital Hand Dynamometer

43100 300 lbs (136 kg) Digital Gauge

#### Digital Pinch Gauge

43101 100 lbs (45 kg) Digital Gauge

#### Inclinometer

43062 Digital Inclinometer Includes carrying case.

43061 **Bubble Inclinometer** 

#### Goniometers

A 43119 Stainless Steel 180° Digit 4" (10 cm)

B 43118 Plastic 180° Digit 4" (10 cm)

# Baseline® Digital Hand Dynamometer

#### Digital Hand Dynamometer

- · Patient will not feel the handle move
- Dynamometer accurately measures grip strength up to 300 lbs (136 kg)
- Handle adjusts to accommodate any hand size
- · Indicator remains at maximum reading until reset
- Includes carrying case
- One-year warranty



# Baseline® Digital Pinch Gauge

#### **Features**

#### **Digital Pinch Gauge**

- Reliably and accurately register up to 100 lbs (45 kg)
- · Measures all pinch tests (tip, key and palmer)
- · Includes carrying case
- One-year warranty



# Baseline® Inclinometer

### Inclinometer

- · Allows total range of motion to be accurately measured
- · Digital model has convenient LED display and carrying case
- Our most popular model, the Bubble Inclinometer is simple and easy to use
- One-year warranty



# Baseline® Goniometers

### **Features**

# Goniometers

- The plastic and stainless steel 180° scale digit goniometers (F and G) require the use of only one hand to measure 110° flexion and 40° hyperextension in 5° increments
- · All goniometers have a linear scale in inches and centimeters
- · One-year warranty



# 3-Piece Baseline® Digital Hand Evaluation Set

### 3-Piece Digital Hand Evaluation Set

- · All the instruments necessary to evaluate strength and range of motion for hand and finger
- · Digital Hand Dynamometer,
- Digital Pinch Gauge (50 lbs/23 kg)
- · Stainless steel finger goniometer
- · Includes carrying case
- · One-year warranty



Ordering Information		
Item#	Description	
3-Piece Bas	eline Digital Hand Evaluation Set	
43102	300 lbs (136 kg) Digital Gauge Dynamometer	
7-Piece Bas	eline Digital Hand Evaluation Set	
43103	300 lbs (136 kg) Digital Gauge Dynamometer 60 lbs (27 kg) Pinch Gauge	
Baseline Or Pinch Gauge	thopedic and Sports Medicine Mechanical e	
43052	30 lbs (14 kg) Dial Gauge (Blue)	
43053	60 lbs (27 kg) Dial Gauge (Red)	

# 7-Piece Baseline® Digital Hand Evaluation Set

#### **Features**

#### 7-Piece Digital Hand Evaluation Set

- · All instruments needed for a full hand evaluation
- Digital Hand Dynamometer,
- 2-point discriminator with 3rd point
- Finger circumference gauge and finger motion gauge
- · Wartenburg pinwheel
- · Mechanical pinch gauge registering up to 60 lbs (27 kg)
- · Stainless steel goniometer
- Digit-O-Meter
- One-year warranty



# Baseline® Orthopedic and Sports Medicine Mechanical Pinch Gauge

# Orthopedic and Sports Medicine Mechanical Pinch Gauge

- · Measures tip, key and palmer pinch strength of healthy hands
- · Results are consistent with the published Mathiowetz studies
- · Indicator remains at the subject's maximum reading until reset
- · Strength readings are in pounds and kilograms up to 60 lbs (27 kg)
- · Lightweight anodized aluminum with wrist strap
- · Includes carrying case
- · One-year warranty



Ordering Information			
Item#	Size	Chest	Bicep
Sully® Should	er Support		
660087-001	X-Small	34"-40"	9"-10½"
660087-002	Small	34"-40"	11"-13"
660087-003	Medium	38"-42"	13"-15"
660087-004	Large	42"-46"	14"-17"
660087-005	X-Large	45"-51"	15"-17"

# Sully® Shoulder Support





#### **Features**

# Sully® Shoulder Support

On the outside, elastic straps attach with hook to the neoprene at any point, in any direction, and with any amount of force you choose. You can functionally stabilize, assist or restrict movement according to the specific needs of each athlete.


On the inside, the perforated, breathable, rubber-like neoprene grips the skin so the support moves with the athlete. It's as if the stabilization straps are attached directly to the body.

The Sully comes complete with a highly adaptable, bifurcated neoprene strap, and two additional 4" neoprene straps. Worn by professional and amateur athletes in all sports! Measure biceps in a non-flexed position.

### **Multiple Applications**



Posterior Instabilities

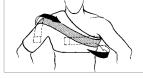

Muscle Strains (pectorial etc.)



Rotator Cuff Deceleration













Anterior Instabilities

Inferior Instabilities

Multi-Directional Instabilities

# Gels and Lotions











### **Features**

### A Ultra-Myossage® Lotion

Formulated with aloe vera, this top quality lotion is ideal for massage or as an ultrasound coupling agent. No menthol, hexachlorophene or polysorbates. Hypoallergenic, water-soluble and greaseless.

### **B** Myossage® Lotion

Soothing, effective alcohol-free, non-greasy and non-staining lotion.

### **C** Utility Straps

A versatile black webbed strap, available in three lengths. Use as a transfer belt, an extension for the EZ-ON traction harness belts, or as shown to assist manual cervical traction.

# **D** Nylatex Wraps

- Ideal for holding cold packs in place
- · Several sizes available
- Reusable

Item#	Description
 A Ultra-M	lyossage® Loti

ion 4262 1 Gallon (3.8 liter) Bottle 4260 8.5 oz. (250 ml) Bottle

B Myossage® Lotion

Ordering Information

1 Gallon (3.8 liter) Bottle 8.5 oz. (250 ml) Bottle (12/case)

C Utility Straps

631374-001 Utility Strap 48" 631325-001 Utility Strap 72" 630483-001 Utility Strap 96"

D Nylatex Assortment (1 pack of 9 rolls)

Nylatex Assortment

D Nylatex Individual (3 rolls per package) \*2.5" x 18" (6 cm x 46 cm) 1264 2.5" x 24" (6 cm x 61 cm) 1204 \*2.5" x 36" (6 cm x 91 cm) \*2.5" x 48" (6 cm x 122 cm) 1208 \*4" x 18" (10 cm x 46 cm) \*4" x 36" (10 cm x 91 cm) \*4" x 48" (10 cm x 122 cm) 4" x 60" (10 cm x 152 cm) 1214 \*6" x 18" (15 cm x 46 cm) \*6" x 36" (15 cm x 91 cm) \*v6" x 48" (15 cm x 122 cm) 6" x 60" (15 cm x 152 cm)

\* Included in Nylatex Assortment Pack (1200)

#### Ordering Information Item# Description DURA-STICK® Premium Self-Adhesive Electrodes 42170 Dura-Stick Premium 1.25" Round 42171 Dura-Stick Premium 2" Round 42172 Dura-Stick Premium 1.5" x 2.5" Oval Dura-Stick Premium 3" x 5" Oval 42173 (20/case) (10 packs of 2) Dura-Stick Premium 2" Square 42174 Dura-Stick Premium 2" x 3.5" Rectangle 42175 42176 Dura-Stick Premium 2" Square Blue Gel Dura-Stick Premium 1.5" x 3.5" 42177 Rectangle Blue Gel **DURA-STICK Self-Adhesive Electrodes** Dura-Stick 2" x 4" Oval 42160 42178 Dura-Stick 2" Square 42179 Dura-Stick 2" x 3.5" Rectangle

Note: all quantities are 10 packs of 4 each, 40 per case, unless otherwise stated

# Dura-Stick® Premium Self-Adhesive Electrodes



- **Dura-Stick® Premium Self-Adhesive Electrodes**
- 10 packs of 4 each, 40/case unless otherwise noted
- Cloth Backing
- · Stainless steel mesh
- · Choice wire and molded connectors
- · Provides optimal current dispersion with low impedance

# Dura-Stick® Self-Adhesive Electrodes



#### **Features**

#### **Dura-Stick® Self-Adhesive Electrodes**

- 10 packs of 4 each, 40/case
- · Carbon silver film layer incorporated into a low profile design
- · Cloth backing for increased "stretchability" and memory
- · Innovative polymer for improved adhesion and durability

# Dura-Stick® Plus Self-Adhesive Electrodes



### **Features**

- Our most economical electrode.
- Permagel™ Hydrogel for superior conformability to body contours
- Blue foam backing
- · Carbon conductive film incorporated into a flexible spun lace backing
- Flexible and sturdy design for durability and ease of use

# Dura-Stick® Disposable Electrodes



# **Features**

# **Dura-Stick® Disposable Electrodes**

- 40/Package
- Multi-day" electrode for extended use
- May be used for 2-5 days
- Solid gel center surrounded by adhesive
- · Beige mesh backing, looks like a bandage
- · Dispose of when finished

Ordering	Ordering Information	
Item#	Description	
Dura-Stick®	Plus Self-Adhesive Electrodes	
42045	Dura-Stick Plus 2.75" x 5" Rectangle	
42159	Dura-Stick Plus 1.5" x 2.5" Oval	
42181	Dura-Stick Plus 2" x 3.5" Rectangle	
42182	Dura-Stick Plus 2" Round	
42183	Dura-Stick Plus 2" Square	
42160	Dura-Stick Plus 2" x 4" Oval	
DURA-STIC	K Disposable Electrodes	
42180	Dura-Stick 2.25" x 2.5" Disposable 40/Package	

Note: all quantities are 10 packs of 4 each, 40 per case, unless otherwise stated

# Ordering Information

Item# Description

# DURA-STICK sEMG Electrode

Dura-Stick sEMG Round White (100 pk/individual)

# Conductive Rubber Electrodes

2852 Black Rubber Carbon 3" (8 cm) Round (2X)

72853 Red Rubber Carbon 3" (8 cm) Round (2X)

# Dura-Stick® sEMG Electrode



#### **Features**

#### **Dura-Stick® sEMG Electrode**

- Round white
- Foam Backing
- 100 pk/individual

# Conductive Rubber Electrodes



### **Features**

#### **Conductive Rubber Electrodes**

- Conductive rubber electrodes
- 3" (8 cm) round
- · Available in black or red

# **Education Manuals**



Manuals to help educate about back pain and daily activites for preventative measures.

# **Features**

#### **Managing Back Pain**

# By H. Duane Saunders, MS PT & Michael S. Melnik, MS, OTR

A photo-filled guide to daily activities for back pain sufferers. Helps your patients in acute pain to safely participate in essential activities such as dressing, child care, kitchen activities, and more. 31 pages.

### **Your Healthy Back**

# By Richard A Pollock, CSP, H. Duane Saunders, MS PT and Michael S. Melnik, MS OTR

The complete guide for back injury prevention, body mechanics, work station design and claims management for supervisors. 24 pages.

#### **For Your Back**

### By H. Duane Saunders, MS PT

The most comprehensive of our manuals. It thoroughly covers back facts, simple anatomy, and causes and treatments for the most common types of back problems. 48 pages.

### **For Your Neck**

### By H. Duane Saunders, MS PT

Discusses common causes of neck injuries and explains preventive strategies. 48 pages.

### **Save Your Back**

# By H. Duane Saunders, MS PT and Michael S Melnik, MS OTR

The photographic guide to patient handling techniques. Essential for health care provider training. 48 pages.

# Ordering Information

#### Item# Description

### Managing Back Pain

630309-001 Managing Back Pain

#### Your Healthy Back

621060-001 Your Healthy Back

#### For Your Back

630308-001 For Your Back

#### For Your Neck

630310-001 For Your Neck

#### Save Your Back

621059-001 Save Your Back





**DJO, LLC** | **A DJO Global Company T** 800.592.7329 **D** 760.727.1280 **F** 800.242.8329

1430 Decision Street | Vista, CA 92081-8553 | U.S.A.

DJOglobal.com/chattanooga

©2015 DJO, LLC 00-2456 Rev.C