

Ordering Information

KleenSpec® Vaginal Specula—Premium 590 Series

- 590XS KleenSpec® 590 Series Premium Disposable Vaginal Specula, X-Small (24/box; 4 boxes/case)
- 59000 KleenSpec® 590 Series Premium Disposable Vaginal Specula, Small (24/box; 4 boxes/case)
- 59001 KleenSpec® 590 Series Premium Disposable Vaginal Specula, Medium (24/box; 4 boxes/case)
- 59004 KleenSpec® 590 Series Premium Disposable Vaginal Specula, Large (18/box; 4 boxes/case)



KleenSpec® Vaginal Specula—Premium 590 Series with Smoke Tube

- 59005 KleenSpec® 590 Series Premium Disposable Vaginal Specula with Smoke Tube, Small (12/box; 4 boxes/case)
- 59006 KleenSpec® 590 Series Premium Disposable Vaginal Specula with Smoke Tube, Medium (12/box; 4 boxes/case)



KleenSpec® Cordless Illumination System

- 79910 Complete KleenSpec® Cordless Illuminator with Charging Station
- 79915 Same as above, for export purposes—specify country and voltage; conforms to IEC 601-1
- 79920 Complete KleenSpec® Cordless Illuminator with Charging Station and Auxiliary Power Cord
- 79925 Same as above, for export purposes—specify country and voltage; conforms to IEC 601-1
- 79900 KleenSpec® Cordless Illuminator
- 73910 Charging Station for Cordless Illuminator; For export purposes: 73912 (European), 73914 (UK), 73916 (Australia)
- 75910 Auxiliary Power Cord Accessory for Cordless Illuminator; For export purposes: 75912 (European), 75914 (UK), 75916 (Australia)
- 59010 Disposable Sheath for Cordless Illumination System (100/box; 5 boxes/case)



In addition, we offer a complete line of vaginal specula and accessories to meet your needs:

KleenSpec® Vaginal Specula—Premium 580 Series w/Sheath

- 58000S KleenSpec® 580 Series Premium Disposable Vaginal Specula with Sheath, Small (18/box; 4 boxes/case)
- 58001S KleenSpec® 580 Series Premium Disposable Vaginal Specula with Sheath, Medium (18/box; 4 boxes/case)
- 58004S KleenSpec® 580 Series Premium Disposable Vaginal Specula with Sheath, Large (12/box; 4 boxes/case)



KleenSpec® Vaginal Specula—Economy Series

- 58600 KleenSpec® 586 Series Economy Disposable Vaginal Specula, Small (25/box; 4 boxes/case)
- 58601 KleenSpec® 586 Series Economy Disposable Vaginal Specula, Medium (25/box; 4 boxes/case)
- 78600 Expendable Illuminator (6/box; 6 boxes/case)—Can be sterilized 50 times before disposal

KleenSpec® Corded Illumination System

- 78810 KleenSpec® Illumination System (with cord)—4.7 V Halogen Light Handle with on/off switch and transformer. 110-130 V, 60 Hz for North America. UL Listed, CSA Certified
- 78815 Same as above, for export purposes—Specify country and voltage; conforms to IEC 601-1
- 78800 Permanent 4.7 V Halogen Light Handle with On/Off Switch
- 73410 Transformer with 5' Cord—110-130 V, 50 Hz. UL Listed, CSA Certified
- 73315 Same as 73410, but for export purposes—Specify country and voltage; conforms to IEC 601-1
- 08800-U 4.6 V to 4.7 V Halogen Replacement Lamp for 78810 Illumination System
- 07800-U 6 V Halogen Replacement Lamp for 78010 Illumination System
- 58010 Disposable Sheath for Corded Illumination System (25/box; 5 boxes/case)



KleenSpec® Storage Systems

- 78820 KleenSpec® Illumination Storage System
- 58014 KleenSpec® Vaginal Speculum Dispenser

Green Series™ Exam Light IV—Another LED innovation from Welch Allyn

- 48810 GS Exam Light IV with Table/Wall Mount
- 48950 Mobile Stand for GS Exam Light IV
- 52640 Disposable Sheath for Exam Light III and GS Exam Light IV (25/box, 5 boxes/case)



For more info, go to www.welchallyn.com/greenseriest


4341 State Street Road, PO Box 220, Skaneateles Falls, NY 13153-0220 USA
(p) 800.535.6663 (f) 315.685.3361 www.welchallyn.com

© 2012 Welch Allyn SM2582 Rev G



WelchAllyn

Advancing Frontline Care™

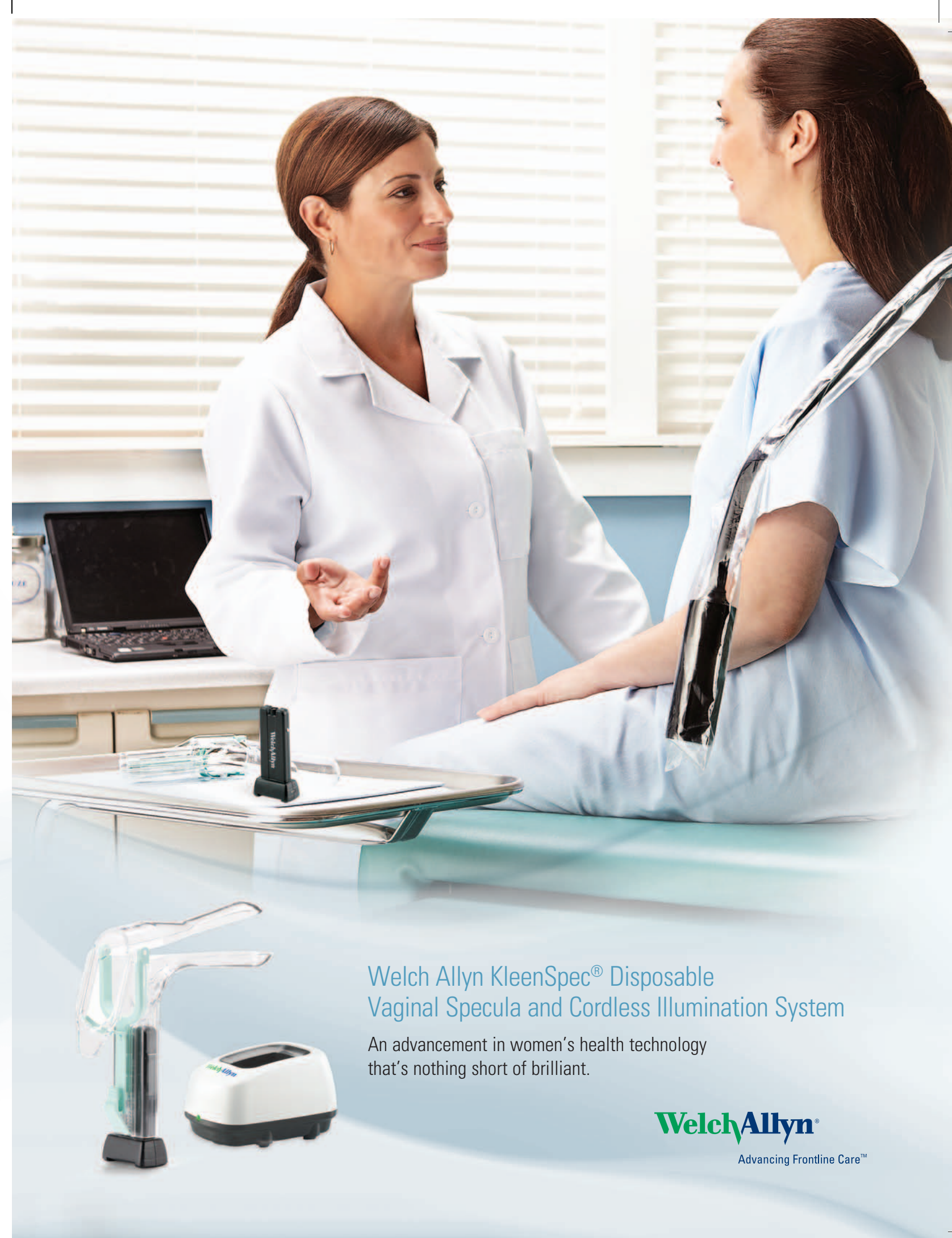

Welch Allyn KleenSpec® Disposable Vaginal Specula and Cordless Illumination System

An advancement in women's health technology that's nothing short of brilliant.

WelchAllyn

Advancing Frontline Care™

Welch Allyn KleenSpec® Disposable Vaginal Speculum and Cordless Illumination System— now in every size needed for all types of patients.



Introducing
the industry's first extra small disposable vaginal speculum

The Welch Allyn 590 Series KleenSpec Disposable Vaginal Speculum

- The new x-small speculum allows you to perform a comfortable exam with virginal/pediatric patients as well as post-menopausal and post-hysterectomy women—no more referrals or uncomfortable exams using metal specula or the wrong size
- Single-use speculum helps to reduce cross-contamination for improved patient safety while eliminating the time and expense of reprocessing metal specs
- Features a wider handle design for improved maneuverability and better balance during exams
- Smooth, molded acrylic designed to prevent pinching or binding of tissue for greater patient comfort
- A strong grip and a yoke that's more flexible than low-cost substitutes

The 590 Series KleenSpec® Vaginal Speculum is compatible with the cordless illumination system (#79910) as well as the corded illumination system (#78810).

Available in x-small (#590XS), small (#59000), medium (#59001) and large (#59004) specula



#59004 Large #59001 Medium #59000 Small #590XS Extra Small



"I had been using the corded version of the Welch Allyn Illumination System since 1980 and was thrilled when the cordless system was introduced. We have been using the Welch Allyn cordless Illumination system with disposable vaginal speculums for all of our patient exams since 2008 and couldn't be happier. This light is brighter than any I have ever seen and shows true tissue color. If you forget to charge the light, it is easy to grab a replacement from another room or use a spare. The light is small, easy and quick to use which truly makes us more efficient. With Manhattan space at a premium, not having floor space taken up by a light source is a great benefit to us. I would recommend this system to any clinicians looking to improve the quality of their care."

Shelley Kolton, MD | Downtown Women OB/Gyn Associates | New York City, NY



The KleenSpec Cordless Illumination System

- LED technology as opposed to halogen, for a whiter, brighter light and true tissue color during an exam
- No bulb replacement necessary for a reduced cost of ownership
- Features a cordless design so you'll never have to worry about cords getting in your way, breaking or needing to be cleaned, speeding the exam and giving you a full range of motion
- The illumination system easily fits into every size speculum
- Newly designed light pipe provides a homogeneous spot and eliminates glare back so all of the light is projected forward
- Uniform spot—no dark or hot spots—for enhanced visualization of the exam area
- Produces much less heat than traditional lamps and has a continuous on time of 80 minutes
- A lithium ion battery provides a recharge time of six hours
- Automatic on/off—insert the illumination system into the vaginal speculum and it automatically illuminates, remove it and the light goes off

