




stretcher series
500
pediatric stretcher/crib

Pedigo's 500 Pediatric Crib Stretcher, with its full function capabilities represents the standard of excellence in hospital cribs. The aesthetically pleasing design makes it perfect for any facility.



Side Rails

Side rails have two-inch bar spacing and a special epoxy finish that is warmer to the touch than chrome.



Total Accessibility

The unique rounded corners reduce the overall size of the 500, yet provide the same amount of useful area on the 26½" x 55¾" mattress.

The two-post support system provides 97% access to the patient by the medical team.



Casters

As a pediatric transport stretcher, the 500 rolls easily on 5" oversized casters. Two of the four casters lock to secure the stretcher while a third steering caster may be engaged for increased control and maneuverability.



A True One-Handed Crib

For complete containment of the child, simply raise the counter balanced side rails to enclose the crib. Because the rails are counter balanced, there will be no sudden drop when released, preventing injury to the child or provider.

Dual-action release handle is childproof and easy to operate, even with a child in one arm.



Crib Security Handle



Squeeze Handle Together



Lift Rail Upward or Downward



STANDARD FEATURES


- 200 lbs. weight capacity rating
- Two 5" locking casters and one steering caster for increased maneuverability
- Trendelenburg/Reverse Trendelenburg
- One handed rail operation
- White powder coat finish
- Unique rounded corners
- Rails with 2" spacing between bars
- Two IV pole receptacles
- Full function fowler
- Litter Top: Solid pan
- Color decals: Blue
- Latex-free

SPECIFICATIONS


Model Number	500 stretcher/crib
Patient Surface	26" x 56¼" (66 x 143 cm)
Overall Length	60" (152 cm)
Overall Width	30¼" (77 cm)
Height	
Fixed Litter Top	33¼" (84 cm)
Top of Crib	59½" (151 cm)
With Security Top	76¾" (195 cm)
Litter Positioning	
Fowler Angles	0°, 20°, 30°
Trendelenburg	0°, 5°, 10°
Fowler Angles	-5°
Side Rails Intermediate Stop	13¾" (35 cm)
Casters	5" (13 cm)
Warranty	1 year parts and labor

OPTIONS

Accessory	Accessory	Accessory
Utility Shelf: Utility shelf 5103001	O₂ Holder: Oxygen tank holder 5114001 Oxygen tank holder 5114001-1 (for larger regulators)	IV Pole: IV Pole 5116001
Security Top Windowed top with six mounting clips 5107001 (will accommodate a mist tent - requires 5104002 for mounting) Clear top with six mounting clips 5117001 (will accommodate a mist tent)	Mattress and Pad: Crib side rail bumper pad (pair) 5109002 Pressure Care Mattress 4"5817002 • Mattress meets fire codes TB-117. Mattress with fire barrier layer 3½" 5860003 • Mattress meets fire codes TB-117 & TB-129. • For additional fire layer barrier protection, contact customer service for a complete listing of fire codes met by this mattress. Mattress 3½" (standard) 5860002 • Phthalate free • compliant with HR 4040.	Color Decals: ■ ■ ■ Red Head/Foot End Labels.... 5118001 Blue Head/Foot End Labels ... 5118002 (standard) Green Head/Foot End Labels. 5118003
Adapters Traction & croupette adapter . 5104002 Traction adapter, pair 5119001 (includes mounting kit) Croupette humidifier adapter . 5108001 Infusion pump adapter 5106001		

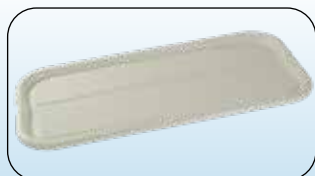

UTILITY SHELF

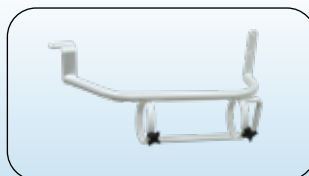

O₂ HOLDERS



CRIB BUMPER PADS