




ProBasics K3 Wheelchair

LIGHTWEIGHT WHEELCHAIR

The **ProBasics K3 Lightweight Wheelchair** offers all the benefits of a standard K3 wheelchair with the added feature of removable flip back height adjustable arms and 2-inch upholstery seat extension. Weighing less than 36 pounds, the ProBasics K3 also comes standard with tool-free adjustable front riggings and easy push-lock wheel locks.



Flip-back armrests are height adjustable

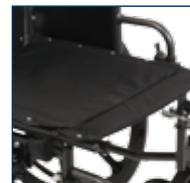

2" upholstery seat extension



Front riggings with heel loop

FEATURES

- Cross-brace, powder-coated silver vein steel frame
- Removable flip back height adjustable arms allow for easy transfer
- 2-inch removable seat extension
- Dual axle provides easy transition of seat height to hemi-level
- Easy-to-clean, comfortable nylon upholstery
- Includes tool-free swing-away footrests with heel loops or elevating legrests with composite footplates
- Elevating legrests are equipped with a pneumatic piston for infinite elevating positions
- Easy push-to-lock wheel locks
- Precision sealed wheel bearings in front and rear ensure long lasting performance and reliability
- Composite handrims with 24" solid urethane tires
- 8" x 1" solid rubber front casters
- HCPCS* Code: K0003

SPECIFICATIONS

Seat Width:	16" 18" 20"
Seat Depth:	16" - 18"
Seat-to-Floor Height ¹ :	Adult 20" Hemi 17"
Back Height:	18"
Back Style:	Fixed
Overall Width w/Arms:	24" 26" 28"
Overall Depth:	31" (without riggings)
Weight Capacity:	300 lb
Product Weight ² :	35 - 38 lb
Carton Dimensions:	33" x 13" x 37"
Shipping Weight:	40 - 46 lb
Warranty:	Limited lifetime on frame One year on all parts

K0195 OR E0990 (elevating legrests only)

Item #	Seat	Arm Style	Leg Style	UPC	UOM
WC31616DS	16" x 16"	Height Adjustable, Flip Back, Desk-Length	Swing-Away Footrests	815067070695	each
WC31816DS	18" x 16"	Height Adjustable, Flip Back, Desk-Length	Swing-Away Footrests	815067070701	each
WC32016DS	20" x 16"	Height Adjustable, Flip Back, Desk-Length	Swing-Away Footrests	815067070718	each
WC31616DE	16" x 16"	Height Adjustable, Flip Back, Desk-Length	Elevating Legrests	815067070725	each
WC31816DE	18" x 16"	Height Adjustable, Flip Back, Desk-Length	Elevating Legrests	815067070732	each
WC32016DE	20" x 16"	Height Adjustable, Flip Back, Desk-Length	Elevating Legrests	815067070749	each

* Providers are responsible for determining appropriate billing codes for Medicare Program claims. This information is included for your convenience and is not intended, nor should it be considered billing advice.

¹NOTE: These heights can vary ¼-inch due to tire wear and manufacturing tolerances.

²NOTE: Weights based on 18" x 16" wheelchair with desk-length arms and no front riggings. Weights will vary depending on how the wheelchair is equipped.