

PRODUCT DATA SHEET

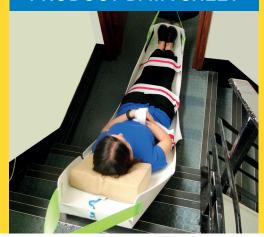

ResQMe

THE RESQME ENABLES CARERS TO RAPIDLY AND SAFELY EVACUATE PEOPLE WITH MOBILITY PROBLEMS IN AN EMERGENCY.

The ResQme is an evacuation board designed to help evacuate non-ambulant people or people who may not be able to walk down stairs in an emergency. It can be used for children to adults 6 foot 6 inches (198cm.) tall. The ResQme can be used in office buildings, theatres, libraries, hospitals, schools and many other facilities.

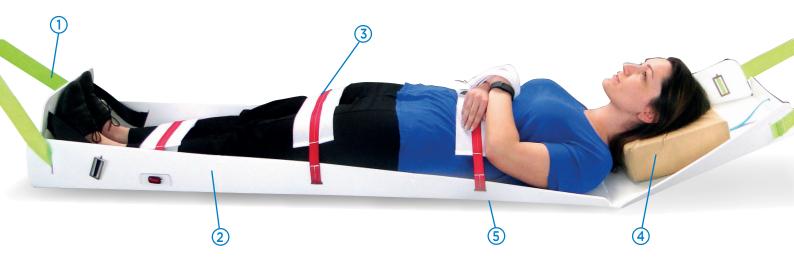

Features

1 Long pulling straps

Located at both head and feet to allow for multiple users

2 Angled side board flaps

Insures the body is kept in place and offering protection when transported around doorways etc.

3 Heavy duty hook and loop

For increased passenger safety

4 Adjustable pillow

For added head protection and allowing to fit persons of varying heights

5 Glide-able plastic board

Corrugated plastic board allows for excellent heat insulation preventing injury, plus removing the necessity to unroll every few months to prevent the plastic keeping its shape

Includes

- 440 lbs Weight capacity
- 2 person operation
- Can be moved easily over virtually any surface
- Narrow width and flexibility
- Wipe clean storage bag
- Can be used down stairs

Benefits

- No installation necessary, ready to use instantly
- No belay system required
- No dangerous ropes to negotiate during evacuation

Dimensions packed

Height	20"
Width	24"
Depth	5"
Weight	4.8lbs

Dimensions unpacked

Height	89"
Width	51"
Depth	3"
Weight Capacity	440lbs

VIDEO LINK: evac-chair.com/evacuation-sleds/resqme/

For over 35 years EVAC+CHAIR® is America's Original and World's Number One Selling Evacuation Chair. Don't Risk Lives, Beware of Poor Imitations and Inferior Substitutions. There is only one EVAC+CHAIR®.

Unsurpassed Quality and Workmanship, Successful in Saving Lives in 9/11, 2 YEAR MANUFACTURERS WARRANTY.

For further information please contact:

Evac+Chair North America LLC

3000 Marcus Ave, Suite #3E6, Lake Success, NY, 11042-1012 **T** 516-502-4240 **W** www.evac-chair.com **E** sales@evac-chair.com

