## **Sports Whirlpools**



Whitehall Manufacturing®

# S Series 110 Gallon Stationary Whirlpool

### S-110-S



Use

### Whitehall S-110-S Stationary Whirlpool is

recommended for treatment and prevention of sports injuries allowing coverage of the extremities, hips and back area. Professional, college and high school athletic trainers find the larger design most suitable when treating today's athletes.

### Tank Capacity is 110 gallons.

**Whirlpool Tank** is 56"L x 24"W x 25"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is seamless, coved (rounded) design for effective cleaning to minimize bacteria buildup. Tank is also provided with an extended stainless steel skirt flush to the floor line to conceal internal drain piping. A 2" Combination Drain and Overflow, and a Filler Spout are provided.

**Turbine Raising and Lowering Device** permits adjustment of desired height and direction of water agitation towards the area of the body to be treated. It is spring balanced and provided with a locking device.

**Turbine Assembly** is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze. **Thermometer** is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

### **GUIDE SPECIFICATIONS**

Provide Whitehall Stationary Whirlpool (specify model number and options) with a 110-gallon capacity, 56"L x 24"W x 25"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Tank shall have a reinforced rim and shall be fitted with a stainless steel skirt flush to the floor line. Bottom of tank shall be a seamless, rounded coved design to minimize bacteria build up. Provide a 2" combination drain and overflow, and a filler spout. Provide a turbine assembly with raising and lowering device which functions as an agitator and permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattle when whirlpool is in operation.

Please visit **www.whitehallmfg.com** for most current specifications.



Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

S-110-S

	AV	
1		
· `		

# Whitehall Manufacturing<sup>®</sup>

Manufacturer of Healthcare and Rehabilitation Products since 1946

115VAC / 60HZ

230VAC / 60HZ

### **MODEL NUMBER AND OPTION SELECTIONS**

Series

S Sports Series

Capacity 🗹 -110 Gallons

Type S Stationary Whirlpool

#### **Turbine Electrical Requirements**

115VAC / 50Hz 230VAC / 50Hz

Country of Use


### For Mixing Valve:

Refer to supplemental MXT15 sheet for details.

# **Sports Whirlpools**

**Suffix Options** (See Whirlpool Options Section for Details)

- -AHC1 Adjustable High Chair with Casters
- -AHC2 Adjustable High Chair with Rubber Tips
- -ARS Arm Rest Support -AS2
  - Adjustable Suspension Seat
- □ -CTS2 Combination Table with Side Seats 🗋 -EG Enviro Glaze Color Powder Coating -Specify Color:\_ (Special Finishes Consult Factory)
- 🗋 -GFI Ground Fault Interrupter
- -PCS Pump Strainer Cover
- Patient Handling System (Chair Type) -PHS2
- 🗋 -RHG **Rubber Hand Grips**
- Separate Drain and Overflow Assembly -SDO
- C -TTS2 Tank Top Seat

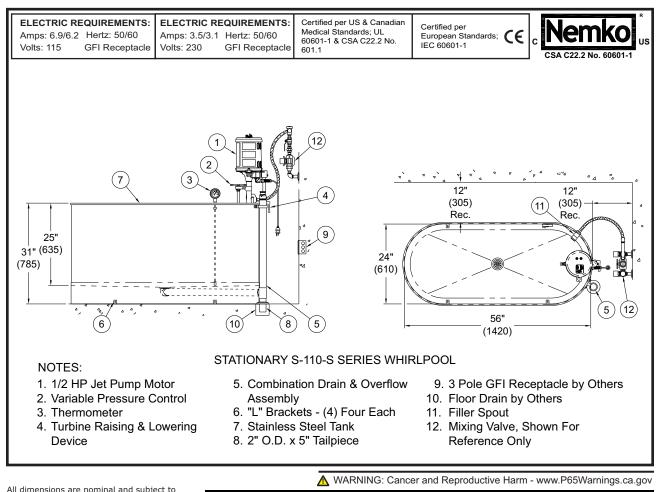

SELECTION SUMMARY & APPROVAL FOR MANUFACTURING

\_Title

manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit www.whitehallmfg.com for most current specifications

Approval for Manufacturing/Signature\_

Model Number & Options\_

Company\_

Contact

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862



Quantity

Date\_