



Cardinal Health™ Protective Apparel

You care for your patients, and we care for you

We are committed to keeping you and your patients safe with our robust portfolio of headwear and footwear that provide the coverage you need, no matter the task.

**CardinalHealth™**

Comfortable coverage



Our headwear and footwear cover all your personal protective equipment (PPE) needs — so you can focus on the task at hand.



Cardinal Health™ Headwear

From surgical caps and hoods to our premium bouffants, our headwear is designed to be comfortable and provide you with necessary coverage at the same time. Cardinal Health™ Headwear is available in a variety of styles, colors and sizes to help you feel most protected.



Cardinal Health™ Footwear

Our durable footwear is constructed of materials ranging from basic spunbond polypropylene to polyethylene-coated spunbond/meltblown/spunbond (SMS) to accommodate a wide range of intended uses. Cardinal Health™ Footwear is available in a variety of styles and sizes to meet your needs.


All Cardinal Health™ Headwear and Footwear is not made with natural rubber latex.

Materials

Spunbond polypropylene

Standard, cost-effective and comfortable


For basic infection control, this non-woven fabric bonds fibers together to form a single layer that is appropriate only for very minimal fluid exposure.

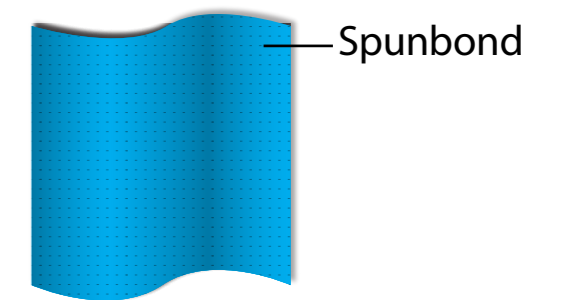

SMS material


Balanced mix of protection and comfort

Strong and breathable, spunbond/meltblown/spunbond (SMS) is a multi-layer fabric composed of inner layers of meltblown polypropylene between outer layers of spunbond polypropylene that is ideal for extended wear. Light/medium weight SMS is appropriate for low amounts of fluid and heavy weight SMS may be appropriate for moderate contact with fluids.

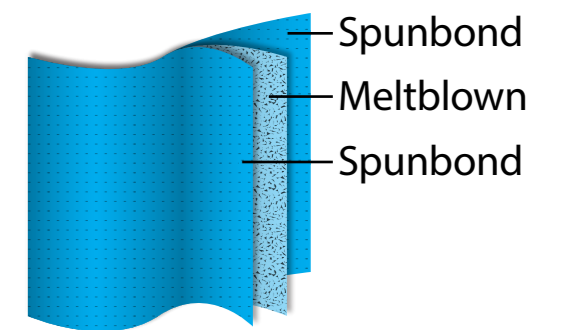

Coated polypropylene (poly-coated)

Soft, spunbond polypropylene is coated with a layer of polyethylene (plastic) film.

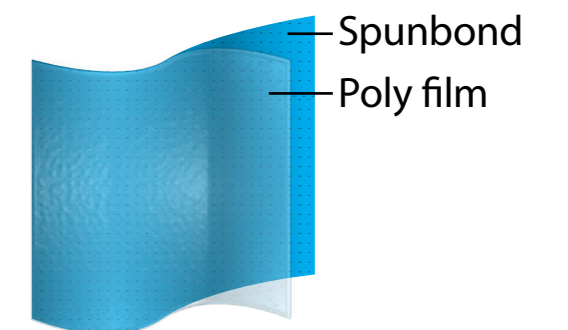


Comfort bouffant caps

Our lightweight, comfortable bouffant caps are made of nonwoven, spunbond polypropylene fabric. The gentle synthetic elastic provides a generous fit and secure coverage without binding.

Cat. No.	Description	Packaging
3224R	Red comfort bouffant, 24 in.	100/bx, 4 bx/cs
3224G	Green comfort bouffant, 24 in.	100/bx, 4 bx/cs
3244	Blue comfort bouffant, 21 in.	100/bx, 10 bx/cs
3274	Blue comfort bouffant, 24 in.	75/bx, 10 bx/cs
3287	Blue comfort bouffant, 27 in.	75/bx, 10 bx/cs
9302	White comfort bouffant, 21 in.	250/cs
9303	White comfort bouffant, 24 in.	250/cs
9304	White comfort bouffant, 27 in.	250/cs
3276	White, celebration print bouffant, 25 in.	75/bx, 10 bx/cs



9303

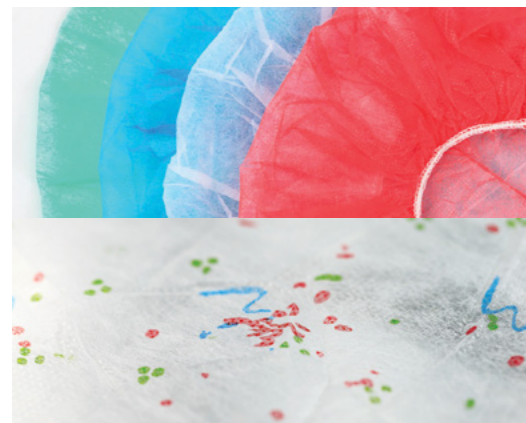

Available in blue, green, red, white and celebration print



Cardinal Health

Premium bouffant caps

Our premium bouffant caps are made of heavyweight spunbond polypropylene. They are designed to provide additional strength while still offering breathability and comfort. Synthetic elastic offers a personalized fit.

Cat. No.	Description	Packaging
B24HD	Blue bouffant, 24 in.	75/bx, 10 bx/cs
3474	Assortment includes smiley faces, bears and tropical print bouffants, 24 in.	75/bx, 10 bx/cs



Available in prints: bears, palm trees and smiley faces



B24HD



3474



Designer scrub cap with ties

A comfortable alternative to bouffant caps, this cap is constructed of lightweight spunbond polypropylene. Designer scrub cap with ties features popular celebration print.

Cat. No.	Description	Packaging
3290	White, celebration print designer scrub cap with ties	100/bx, 6 bx/cs



3290

Surgical cap

Easy-tie surgical cap is designed for particulate control. Our cap is easy to don and features adjustable ties. Lightweight, nonwoven material is comfortable and absorbent.

Cat. No.	Description	Packaging
4359	Blue surgical cap with ties	100/bx, 6 bx/cs
4359SMS	Blue SMS surgical cap with ties	100/bx, 5 bx/cs



4359SMS

Extended surgical cap

Our extended surgical cap benefits your patients by reducing the risk of exposure to potentially pathogenic organisms that live on the hair, skin, ears and facial hair of perioperative team members. This single-use cap has a breathable top, fold-to-fit personalized coverage and two securement options to meet your comfort needs and preferences.

Cat. No.	Description	Packaging
4360	X-Large surgical cap with ties, extended coverage	600/cs
4369	Large surgical cap with elastic back, extended coverage	600/cs
4370	X-Large surgical cap with elastic back, extended coverage	600/cs



Comfort hood

Our cool and comfortable hood conforms to the head and neck, providing maximum coverage for hair and sideburns. Adjustable ties at the neck offer a secure and comfortable fit.

Cat. No.	Description	Packaging
4380C	Blue comfort hood	75/bx, 6 bx/cs



4380C

Surgical pullover hood

This surgical hood provides complete head coverage. Blue hood has gentle elastic at neck for a comfortable fit.

Cat. No.	Description	Packaging
4384	Blue pullover hood	100/cs



4384

Surgical hood with ties

Lightweight and comfortable, our hood provides complete head, neck and shoulder coverage. Adjustable back ties allow for a personalized fit.

Cat. No.	Description	Packaging
4381	Blue surgical hood with ties	100/cs



4381

Surgical beard cover

White beard cover provides comfortable coverage. Elastic band fits comfortably around head.

Cat. No.	Description	Packaging
9216	White beard cover, full coverage	250/cs



9216



Spunbond polypropylene shoe covers

Our lightweight, comfortable shoe covers are made of nonwoven, spunbond polypropylene fabric. Covers are available with or without skid-resistant tread.

Cat. No.	Description	Size	Packaging
2850	Basic shoe cover	Universal	100/bx, 4 bx/cs
2852	Skid-resistant shoe cover	Universal	100/bx, 4 bx/cs
2854	Skid-resistant shoe cover	X-Large	100/bx, 4 bx/cs



SMS shoe covers

Cardinal Health™ Dura-Fit™ and Dura-Pak™ shoe covers are made of SMS fabric. Covers are available with or without skid-resistant tread.

Cat. No.	Description	Size	Packaging
4850	Dura-Fit™ Shoe Cover	Universal	100/bx, 4 bx/cs
4852	Dura-Fit™ Skid-Resistant Shoe Cover	Universal	100/bx, 4 bx/cs
4854	Dura-Fit™ Skid-Resistant Shoe Cover	X-Large	100/bx, 4 bx/cs
4855	Dura-Fit™ Skid-Resistant Shoe Cover	XX-Large	100/bx, 4 bx/cs
4859SS	Dura-Fit™ Skid-Resistant Shoe Cover	Sport	100/bx, 4 bx/cs
4872A	Dura-Pak™ Skid-Resistant Shoe Cover	Universal	400/cs
4874	Dura-Pak™ Skid-Resistant Shoe Cover	X-Large	400/cs
4875	Dura-Pak™ Skid-Resistant Shoe Cover	XX-Large	400/cs



Heavy-duty SMS shoe covers

Our extra-strong SMS shoe covers provide increased strength. Covers are available with or without skid-resistant tread.

Cat. No.	Description	Size	Packaging
5850	Shoe cover	Universal	100/bx, 4 bx/cs
5852	Skid-resistant shoe cover	Universal	100/bx, 4 bx/cs
5854	Skid-resistant shoe cover	X-Large	100/bx, 4 bx/cs
5859	Skid-resistant shoe cover	XX-Large	100/bx, 4 bx/cs



Coated polypropylene (poly-coated) shoe covers

For use in fluid-intensive conditions, Cardinal Health™ Critical Zone™ and FullGuard™ Shoe Covers are made of durable polypropylene fabric that meets ASTM F1670 and ASTM F1671. Patented separate elastic band at the ankle minimizes the chance of fluid strike-through by eliminating pinholes caused by stitched-in elastic. Bound elastic seams at the foot reinforce fluid protection and adjust for a customized fit. Foam strips provide anti-skid protection and assist in donning orientation. Critical Zone™ Shoe Covers offer protection in the shoe area, while FullGuard™ Shoe Covers provide fluid-intensive protection throughout.

Cat. No.	Description	Size	Packaging
8357	Ankle-top shoe cover	Universal	30/bx, 4 bx/cs
8358A	Ankle-top shoe cover	X-Large	30/bx, 4 bx/cs
8451	Economy Zone™ High-Top Shoe Cover, glue-based tread	Universal	20/bx, 4 bx/cs
8452	Economy Zone™ High-Top Shoe Cover, glue-based tread	X-Large	20/bx, 4 bx/cs
8453	Critical Zone™ High-Top Shoe Cover	Universal	20/bx, 4 bx/cs
8454	Critical Zone™ High-Top Shoe Cover	X-Large	20/bx, 4 bx/cs
8457	FullGuard™ High-Top Shoe Cover	Universal	20/bx, 4 bx/cs
8458	FullGuard™ High-Top Shoe Cover	X-Large	20/bx, 4 bx/cs



Your headwear and footwear solution for staff safety

For more information or to order, please contact your Cardinal Health sales representative, call **877.CARDINAL** or visit **cardinalhealth.com**

SmartGown
CardinalHealth

