

See the Difference. Feel the Quality.

Uni-Patch™ R Series (Re-Ply) Stimulating Electrodes

Your patients rely on you for care — you can rely on us to help with our Uni-Patch™ R series (re-ply) stimulating electrodes. Designed for durability and comfort, they feature:

- A long-lasting 0.025-inch thick hydrogel for reliable adhesion that conforms well over body contours
- A quick, moderately aggressive tack for firm yet gentle electrode-to-skin contact for multiple applications
- Flexible, durable, and highly conductive pigtail wire strands designed for multiple uses
- Highly conductive carbon film for even current distribution
- Durable tricot covers to help extend electrode life

The soft mold flex connector with secure lead wire adhesive technology features:

- A socket that easily adapts to variations in pin size and provides a secure connection
- Wire adhesive designed to anchor the lead wire into the pad and diminish separation





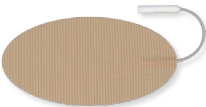



When you're done using the electrodes, you can return them to either side of the double-coated, tinted storage liner.



R Series (Re-Ply)
Stimulating Electrodes

ORDERING INFORMATION

Uni-Patch™ R Series (Re-Ply) Stimulating Electrodes

CODE	MODEL	SIZE	SHAPE	QTY
 EP84770	650	1.5 x 2 in.	Oval	4/pack
 EP84820	652	2 in.	Round	4/pack
 EP84850	653	2.75 in.	Round	4/pack
 EP84910	654	2 x 2 in.	Square	4/pack
 EP84772	657	2 x 4 in.	Oval	4/pack
 EP85040	658	2 x 4 in.	Rectangle	4/pack
Snap connection				
 ES82180	655	2 x 2 in.	Square	4/pack
 ES82200	659	2 x 4 in.	Rectangle	4/pack

10–14 applications

For more information, please contact your Cardinal Health sales representative at OEMSpecialty@cardinalhealth.com or visit cardinalhealth.com/pain-therapy-recovery.

© 2018 Cardinal Health. All Rights Reserved. CARDINAL HEALTH, the Cardinal Health LOGO and UNI-PATCH™* are trademarks of Cardinal Health and may be registered in the US and/or in other countries. All other marks are the property of their respective owners. 02/2018 - H9830 - WF#755854

15 Hampshire Street
Mansfield, MA 02048

800.962.9888
508.261.8000

