

## **Hydrolift**

# HL Series Hydrolift

HL-10, HL-15-M, HL-22-M

(Model HYL-15-M shown in photo includes HL-15-M with P-)



Product photo/image above may show some available options

The Whitehall HL Series Hydrolift is recommended for the raising and lowering of a single Whitehall Hydrotherapy Whirlpool to treat lower and upper extremities. The Hydrolift is fully adjustable, raising the whirlpool rim to any height up to 32", 36" or 38" depending on the model being specified. In normal operating conditions, the Hydrolift will remain as a stationary lift. However, both Hydrolift and whirlpool may be moved by lowering the whirlpool to the floor and using the whirlpool casters.

New or existing Whitehall Hydrotherapy Whirlpool models P-10-M, P-15-M, E-15-M, P-22-M, E-22-M can be used with the HL series Hydrolift. Due to variable dimensions, indicate if whirlpool suffix option -SDP (Separate Drain Pump) is specified with whirlpool.

Hydrolift Frame, Clamp and Push Bar are constructed of heavy gauge, type 304 stainless steel, polished to a satin finish. The base of the frame sits underneath the whirlpool to prevent any possible obstacle. The frame is of a telescoping design and includes a tamper resistant clamp to secure the Hydrolift to the whirlpool with a push bar to assist when moving.

**Crank Handle** to raise and lower whirlpool is easily adjustable and includes an integral brake lever to lock the handle and whirlpool in place. A tension adjuster is provided at the handle to regulate any changes in weight being lifted.

#### **GUIDE SPECIFICATIONS**

Provide a Whitehall Hydrolift (specify model number and options) for use with new or existing Whitehall Hydrotherapy Whirlpool (specify model number). Hydrolift is to raise whirlpool rim above the floor (specify height dimension based on Hydrolift model being specified). Frame is to include a vandal resistant clamp for attaching the whirlpool and push bar to assist when moving the whirlpool. Frame, clamp and push bar are to be steel, with a satin finish. Crank handle is to include an integral brake lever to lock the handle in place and tension adjuster.

Last Update: 2/14/08



# Hydrolift

#### **MODEL NUMBER AND OPTION SELECTIONS**

#### Series

O -HL Hydrolift

#### Capacity

- O -10 Gallons
- O -15 Gallons
- O -22 Gallons

#### **Type**

O -M Mobile

**Suffix Options** (See Whirlpool Option Section for Details)


O -EG Enviro Glaze Color Powder Coating -Specify Color:

(Special Finish Consult Factory)


O -SDP Separate Drain Pump


*	=	Add	2"	on	-SDP	Separate	Drain	Pump.
---	---	-----	----	----	------	----------	-------	-------

MODEL	DIMENSION								
MODEL	Α	В	C *	D *	E	F			
III 10 M	38"	53"	17"	32"	19-1/2"	19"			
HL-10-M	(965)	(1345)	(430)	(815)	(495)	(485)			
111 15 M	38"	53"	20"	35"	19-1/2"	21"			
HL-15-M	(965)	(1345)	(510)	(890)	(495)	(535)			
ш 22 м	41"	56"	23"	38"	21-1/2"	24"			
HL-22-M	(1040)	(1420)	(585)	(965)	(545)	(610)			

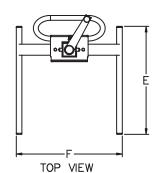

### NOTES:

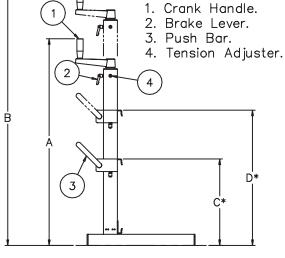

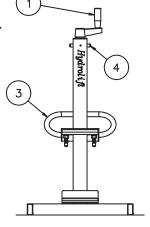

SIDE VIEW

FRONT VIEW

### HL SERIES HYDROLIFT

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility for use of void or suspended data. © 2008 Whitehall Manufacturing, A Division of Acorn Engineering Company. Please visit www.whitehallmfg.com for most current specifications.

#### SELECTION SUMMARY & APPROVAL FOR MANUFACTURING

Model Number & Options\_ Quantity Company\_ Date Title Contact Approval for Manufacturing/Signature

Whitehall Mfg. > P.O. Box 3527 > City of Industry, CA 91744 > (800) 782-7706 > (626) 968-6681 > FAX (626) 855-4862