

The acapella<sup>®</sup> product family.  
Giving patients a breadth of fresh air.




# Greater product breadth to better treat your patients.

The acapella® vibratory positive expiratory pressure (PEP) products are especially helpful in achieving the following therapeutic goals of bronchial hygiene therapy:

- Aiding in mobilization of retained secretions (in cystic fibrosis and chronic bronchitis patients)<sup>1</sup>
- Preventing or reversing atelectasis<sup>1</sup>
- Reducing air trapping in patients with asthma and COPD<sup>1</sup>
- Optimizing delivery of bronchodilators in patients receiving bronchial hygiene therapy<sup>1</sup>

acapella® products offer unique and convenient features that allow time savings for therapists and patients alike. Your patients won't be the only ones breathing easier!

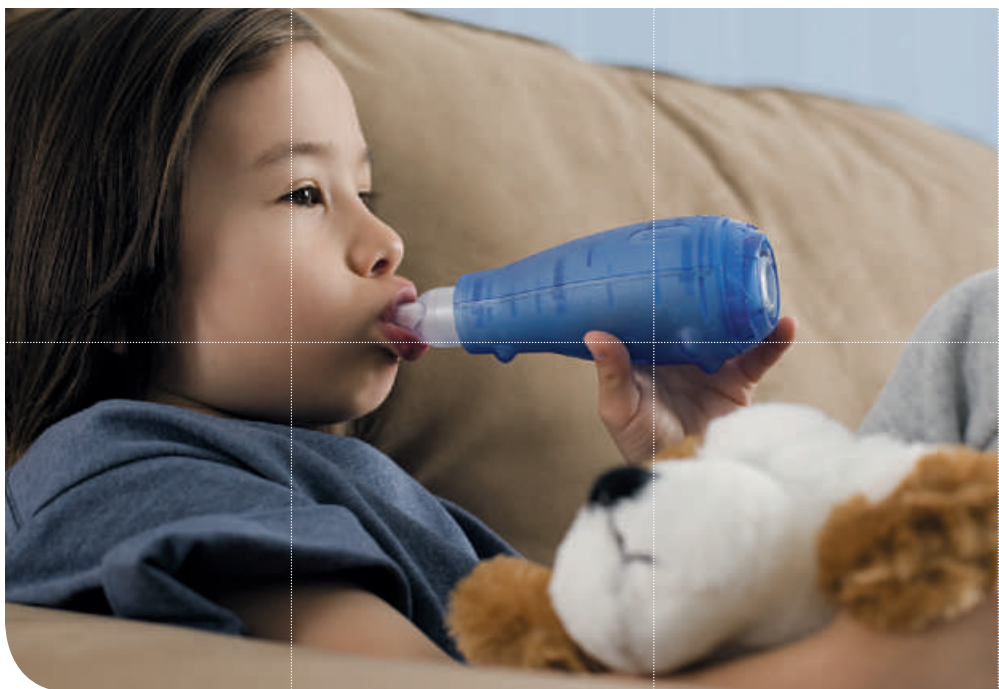


Key Features and Benefits	acapella® duet	acapella® choice	acapella® DH	acapella® DM
Functional in any position – Trendelenberg, standing, or sitting	•	•	•	•
Adjustable frequency and flow resistance settings	•	•	•	•
Treatment can be performed in less than half the time of a conventional CPT session <sup>2</sup>	•	•	•	•
Distal 22mm O.D. fitting allows nebulizer connection via tee adapter		•	•	•
Convenient built-in nebulizer port	•			
Streamlined body design offers easy grip	•			
Tethered cap helps reduce risk of contamination when in resting position	•			
Easily disassembled for heat disinfection via boiling, autoclaving, or placing in the dishwasher (top rack only)	•	•		
≥ 10 L/min expiratory flow requirement	•	•		
≥ 15 L/min expiratory flow requirement			•	
< 15 L/min expiratory flow requirement				•

<sup>1</sup>AARC Clinical Practice Guidelines, Use of Positive Airway Pressure Adjuncts to Bronchial Hygiene Therapy, published in Respiratory Care 1993;38(5):516–521.


<sup>2</sup>Mahlmeister MJ et al. "Positive-expiratory-pressure mask therapy: Theoretical and Practical Considerations and a Review of the Literature," Respiratory Care, 1991.



## acapella® product family features

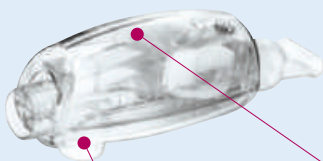

Standard-sized port will accommodate most medication nebulizers  
**(acapella® duet only)**

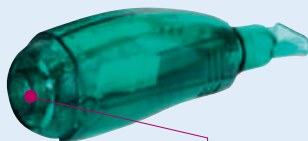

Clear coloring aids in visual recognition of cleanliness  
**(acapella® duet only)**

One-way inspiratory valve allows inhalation without removal from the patient's mouth  
**(all acapella® products)**



Expiratory resistance dial for customized frequency adjustment  
**(all acapella® products)**



Distal fitting allows easy fit with 22mm I.D. connections  
**(all acapella® products)**

Proximal 22mm O.D. connection allows use with a mouthpiece or mask  
**(all acapella® products)**

## The acapella® vibratory PEP product family



acapella® duet




acapella® choice



acapella® DH



acapella® DM

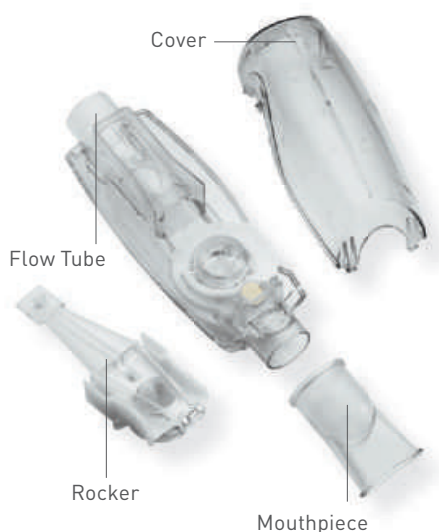

acapella® duet and acapella® choice vibratory PEP products are easily disassembled and reassembled for effortless cleaning.

For more information about the acapella® family of products, visit [www.smiths-medical.com](http://www.smiths-medical.com).

PART #	PRODUCT DESCRIPTION	UNITS/CASE
<b>acapella® duet vibratory PEP products</b>		
27-9000	acapella® duet kit, includes acapella® duet, mouthpiece, Portex® SVN, oxygen tubing, collapsible flex tubing	10
27-9001	acapella® duet and mouthpiece only	10
<b>acapella® choice vibratory PEP products</b>		
27-7000	acapella® choice with mouthpiece	10
<b>acapella® DH (green) vibratory PEP products</b>		
21-1530	acapella® DH with mouthpiece	10
21-3530	acapella® DH with pediatric mask	1
21-5530	acapella® DH with medium mask	1
21-7530	acapella® DH with large mask	1
<b>acapella® DM (blue) vibratory PEP products</b>		
21-1015	acapella® DM with mouthpiece	10
21-3015	acapella® DM with pediatric mask	1
21-5015	acapella® DM with medium mask	1
21-7015	acapella® DM with large mask	1
<b>acapella® accessories</b>		
20-0022	TheraPEP® pressure indicator, tubing & pressure port	10
27-0050	Replacement mouthpiece; fits all acapella® family designs	50

Note: Product is for single patient use only.

PRODUCT(S) DESCRIBED MAY NOT BE LICENSED OR AVAILABLE FOR SALE IN CANADA AND OTHER COUNTRIES

**Smiths Medical International Ltd.**  
 Ashford, UK TN25 4BF  
 Phone: +44 (0) 845 8500445  
[www.smiths-medical.com](http://www.smiths-medical.com)

Find your local contact information at: [www.smiths-medical.com/customer-support](http://www.smiths-medical.com/customer-support)

Smiths Medical is part of the global technology business Smiths Group plc. Please see the Instructions for Use/Operator's Manual for a complete listing of the indications, contraindications, warnings and precautions. Portex, acapella and the Smiths Medical and Portex design marks are registered trademarks of Smiths Medical. The symbol © indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. All other names and marks mentioned are the trademarks or service marks of their respective owners. ©2013 Smiths Medical. All rights reserved. RE194202GB-082013

CE Rx  
0473 ONLY

smiths medical

MROCCA-0015