


Discover

Simplified, Multi-Drug Syringe Delivery in the Operating Room

InfusO.R. Pump with SMART LABELS

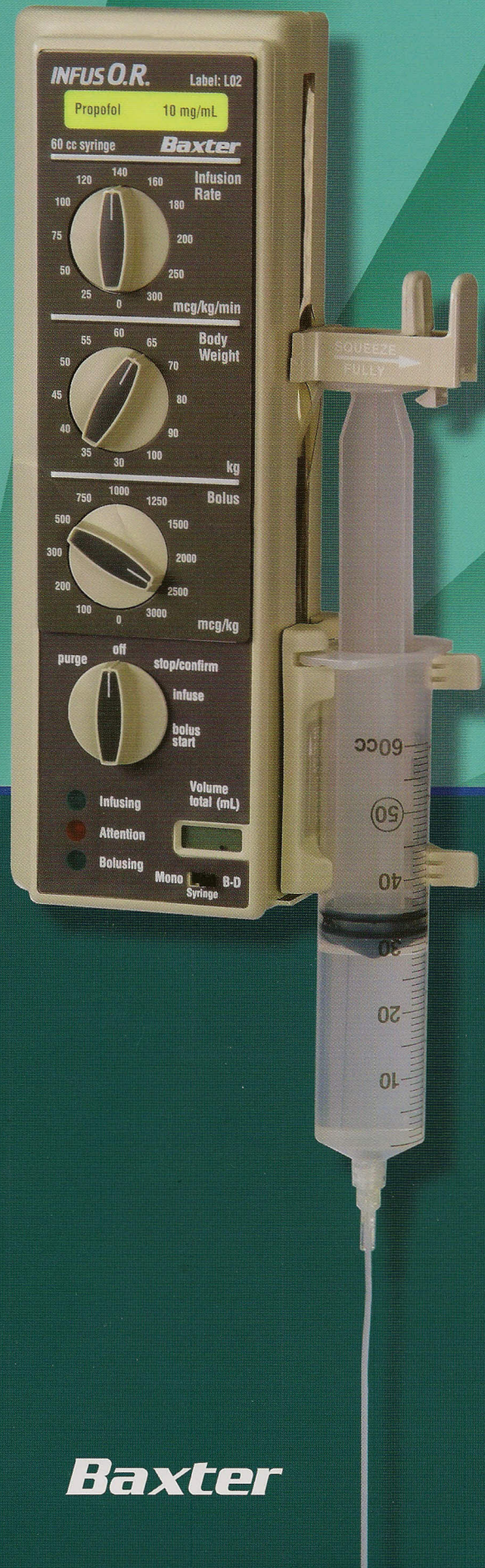

Put Ease-of-Use at Your Fingertips

- Drug specific SMART LABELS* help minimize the potential for medication errors in dose programming
- Programmed settings may be adjusted easily at any time during infusion
- SMART LABELS* are specifically designed to work with either a 20 cc syringe or 60 cc plastic syringe

Bringing Infusion Versatility to the Operating Room

- Designed for use in the OR and for outpatient surgical settings
- Specific SMART LABELS* support Pediatric and Monitored Anesthesia Care infusion therapies

* Each SMART LABEL converts the InfusO.R. Pump into a specific delivery system for the indicated drug with appropriate rate selections and delivery units. There are also SMART LABELS available for General and General-Bolus infusion applications.

Baxter

InfusO.R. Syringe Pump

SPECIFICATIONS

Product Code

2L3100

Accuracy

Infusion: +/- 3% linear rate
Bolus: +/- 3% linear displacement

Battery Life

150 hours typical use and 4 hours at low battery

Batteries

Four C Size alkaline batteries (ANSI/NEDA 14A)

Syringes

Kendall Monoject® or B-D® 20 cc plastic
Kendall Monoject® or B-D® 60 cc plastic

Flow Rates

0-600 mL/hr, depending on Smart Label

Flow Profile

Infusion: pulsed intermittently
Bolus: continuous

Occlusion Force (Pressure)

8 ± 1 lbs

Maximum Occlusion Pressure

11 psi with 60 cc syringe
20 psi with 20 cc syringe

Occlusion Detection Time*

(Inversely related to flow rate)
120 seconds at 36 mL/hr

Volume Stored on Occlusion*

1.1 mL (approximate)

Back Pressure Effect on Accuracy

None to occlusion pressure

Size

9.2 x 4.5 x 2.0 inches (238 x 113 x 50 mm)


Weight

2 lbs (1 kg) with batteries

* Using a microbore anesthesia set and a 60 cc syringe

Rx Only. For safe and proper use of this device, please refer to the complete instructions in the user's manual.

ADMINISTRATION SETS

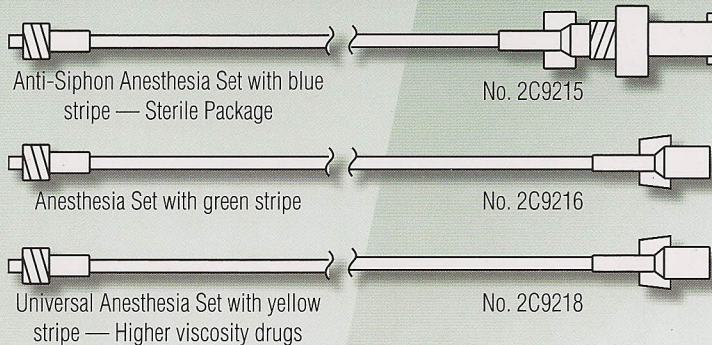

SMART LABELS

Medical Description	Medical Concentration	Catalog Number	20 cc Syringe	60 cc Syringe
Cisatracurium	2 mg/mL	2L3260		✓
Dopamine	4000 mcg/mL	2L3230		✓
Fentanyl	50 mcg/mL	2L3243	✓	
General mL/hr	Not Applicable	2L3241		✓
General mL/hr-Bolus	Not Applicable	2L3225		✓
Mivacurium	2 mg/mL	2L3248		✓
Propofol	10 mg/mL	2L3242		✓
Propofol-Pediatric Dose	10 mg/mL	2L3259		✓
Remifentanil-General Anesthesia	50 mcg/mL	2L3269		✓
Remifentanil-Monitored Anesthesia Care	25 mcg/mL	2L3268		✓
Sufentanil	50 mcg/mL	2L3235	✓	
Sufentanil	5 mcg/mL	2L3249		✓
Sufentanil-Pediatric Dose	10 mcg/mL	2L3240	✓	
Vecuronium	1 mg/mL	2L3227	✓	

**For more information contact your
Baxter Sales Representative,
or call 1-800-933-0303.**



MEDICATION DELIVERY

Baxter, InfusorO.R. and Committed to a Safer Healthcare Environment are trademarks of Baxter International Inc.

B-D and Monoject are registered trademarks of Becton-Dickinson Co. and Tyco Healthcare Group Lp, Inc.

Baxter Healthcare Corporation, Medication Delivery, Route 120 and Wilson Road, Round Lake, IL 60073 www.baxter.com 420255 5K 10/08

Baxter