


Wall-mounted stadiometer with display on the headpiece

seca 264

- + Wall mounting: Absolutely stable and wobble-free height measurement
- + Digital height measurement: comfortable and precise
- + Precise head alignment by the seca Frankfurt line positioner
- + Headpiece with a backlit digital display for easy reading
- + Non-slip rubber mat with heel positioner

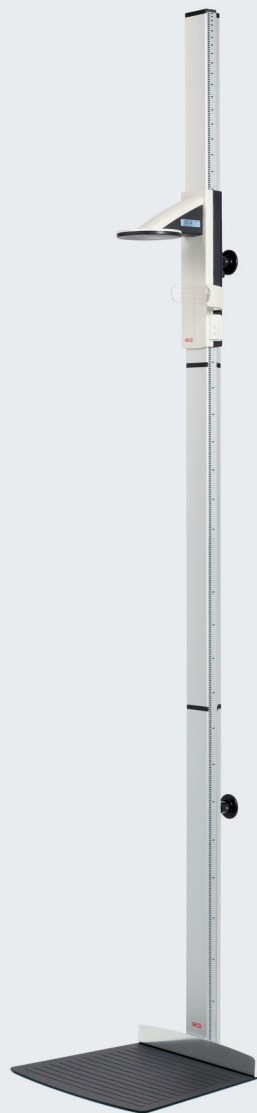

seca 264

Wall-mounted stadiometer with display on the headpiece



Highly stable wall mounting for the most precise height measurement

The customized column made of custom-made aluminum profiles is firmly connected to the wall with two brackets. This gives the seca 264 stadiometer an extremely high level of stability, making measuring height particularly easy and comfortable. The digital height caliper has the seca Frankfurt Line, which assists to achieve ultra-precise head alignment.

Measuring range	11" – 7 ft 2" / 30 – 220 cm
Graduation	1/8" / 1 mm
Functions	HOLD, cm/inch switch-over, user-defined zero setting, automatic switch-off, SEND



Convenient, fast and precise height measurement

This is what the perfect height measurement with the seca 264 stadiometer looks like: Align the patient's feet with the heel positioner and guide the head slide with brake button onto the crown. The head is perfectly aligned horizontally thanks to the Frankfurt Line positioner. The height can be read directly from the display in the head slide.



Ultra-stable wall stadiometer for perfect height measurement

Mounting the seca 264 stadiometer on the wall provides enormous, firm stability. Together with the heel positioner and the seca Frankfurt line positioner, the patient can be precisely positioned and the height can be determined quickly and precisely. The measuring result is read on the display of the headpiece.

We reserve the right to make modifications.

2011264-04022

seca north america | west
13601 Benson Avenue
Chino, CA 91710
USA

seca north america | east
498 Wando Park Boulevard
Mount Pleasant, SC 29464
USA

T +1 800 542 7322
F +1 888 705 7397
E info.us@seca.com

Precision for health
seca.com